



User Manual
External Semi-automatic Defibrillator
for public access

Rev.	4.2
------	-----



QUICK USE GUIDE



©by A.M.I Italia s.r.l.

These instructions for use cannot, without our consent, be completely or partially reproduced, transmitted, stored electronically or translated into another foreign language or computer language. Infractions against this prohibition not only violate our copyright, but also reduce our ability to provide accurate and up-to-date information to the user and the operator of the device.

A.M.I Italia S.r.l.

Via Cupa Reginella, 15A - 80010 Quarto (NA) Italy

Tel. +39 081 806 34 75 Fax +39 081 876 47 69

e-mail info@amiitalia.com

<http://www.amiitalia.com>

Printed in Italy

Summary

1	Introduction.....	8
1.1	Preface.....	8
1.2	Use in accordance with provisions.....	8
1.3	Guarantee.....	8
1.4	Exclusion of liability.....	8
1.5	Indications.....	8
1.6	Counter indications.....	8
1.7	Version information.....	9
1.8	Symbols in the manual.....	9
1.9	Manufacturer contacts.....	9
2	Safety instructions.....	10
2.1	Indications of DANGER.....	10
2.2	Indications of CAUTION.....	10
2.3	Indications of DISPOSAL.....	11
3	Description of the device.....	12
3.1	Device Information.....	12
3.2	Service activation procedure.....	14
3.3	Classifications.....	14
4	Description of device details.....	15
4.1	General structure of device.....	15
4.2	Keys, icons and indicators.....	16
4.3	Status display.....	16
4.4	Standard and optional accessories of the device.....	17
5	Parts and accessories of the Geo Saver.....	19
5.1	AED Batteries.....	19
5.1.1	Non rechargeable Li-SOCl ₂ battery (SAV-C1032).....	19
5.1.2	Suggestions for a proper maintenance of battery SAV-C1032.....	Errore. Il segnalibro non è definito.
5.1.3	Rechargeable Li-Ion battery (SAV-C1033).....	20
5.1.4	Suggestions for a proper maintenance of battery SAV-C1033.....	20
5.1.5	Inserting and removing the batteries.....	21
5.1.6	Recharging station for rechargeable batteries.....	22
5.1.7	Battery charger structure.....	22

5.1.8	Recharge procedure	23
5.2	Geoloc Batteries	24
5.2.1	Non-rechargeable Li-SOCl ₂ battery (SAV-C1038)	24
5.2.2	Rechargeable Li ion battery (SAV-C1039).....	25
5.2.3	Inserting and removing the batteries.....	26
5.2.4	Recharge procedure	27
5.3	Defibrillation PADS.....	28
5.3.1	Defibrillation PADS for Adults SAV-C0846	28
5.3.2	PADS for Children SAV-C0016.....	28
5.3.3	Positioning of defibrillation PADS.....	29
5.4	Memory Card and USB	30
6	Auto test.....	31
6.1	Display and control LED.....	31
6.2	ACTIVATION Test	32
6.3	AUTOMATIC Test.....	33
6.4	POWER ON Test.....	33
7	Information button.....	34
8	Defibrillation.....	35
8.1	Switching on the <i>Geo Saver</i>	35
8.2	Preparation of the patient.....	36
8.3	Positioning the PADS	36
8.4	Cardiac rhythm analysis.....	37
8.5	Shockable Rhythm	38
8.6	Change of rhythm	39
8.7	Non-shockable rhythm	39
8.8	CPR.....	39
9	Recording, viewing and archiving the data	41
9.1	Data recording	41
9.2	Archiving data on PC.....	41
10	Web Management.....	42
10.1	Main features	43
10.1.1	Changing the language	43
10.1.2	User profile and password.....	43
10.1.3	Device management.....	44
10.1.4	ECG streaming	45

10.1.5	Auto tracking	46
10.1.6	Battery replacement.....	46
10.1.7	Firmware update	47
10.2	Dashboard	48
10.3	Devices.....	49
10.4	Geo Saver.....	50
10.5	Saver Log.....	51
10.6	Geoloc Log	52
10.7	Streaming.....	53
10.8	List of saved ECG Streaming	54
10.9	Tracks.....	55
11	Maintenance.....	56
11.1	After each use.....	56
11.2	Ordinary Maintenance	56
11.3	Cleaning	57
11.4	Preservation.....	57
11.5	Troubleshooting Guide	58
12	Technical Specifications.....	59
12.1	Physical characteristics	59
12.2	Environmental requirements	59
12.3	Reference regulations	59
12.4	Alarms Table	60
12.5	Controls and indicators.....	60
12.6	Data memory	60
12.7	Defibrillator.....	61
12.8	Efficiency of delivered energy	62
12.9	Patient analysis system.....	64
12.10	ECG Analysis Function	64
12.11	Non rechargeable battery	64
12.12	Rechargeable battery	65
12.13	Internal back-up battery.....	65
12.14	Rechargeable battery charger	65
12.15	Defibrillation PADS.....	65
12.16	Timing of Shock cycles	66
12.17	Geoloc Module	66

12.18	Geoloc Module Battery Type.....	66
13	Compliance with electromagnetic emission standards.....	67
13.1	Guidelines and manufacturer's declaration - Electromagnetic emissions.....	67
13.2	Guidelines and manufacturer's declaration – Electromagnetic immunity.....	67
13.3	Recommended separation distance between portable and mobile RF communication equipment and Geo Saver device	69
14	Simbology	70
15	Certifications.....	71
15.1	EC Certificate	71
15.2	IMQ Brand	73
16	Geo Saver Series Defibrillator Warranty	74
17	Product registration.....	75

1 Introduction

1.1 Preface

Thank you for having chosen the defibrillator of A.M.I Italia S.r.l. model **Geo Saver**.

So that you can correctly use the device it is necessary, before usage, to carefully read this user manual. The User Manual of **Geo Saver** contains the instructions for its use in compliance with its function and purpose. For a function free of error and to achieve the right benefits, it's fundamental to respect the prescriptions indicated in this user manual, to guarantee the safety of the patient, of the rescuer and of any third parties. This manual is an integral part of the defibrillator and must always be kept together with the device, so that it can be easily accessible if necessary.

1.2 Use in accordance with provisions

The device **Geo Saver** can be used exclusively if the conditions indicated in the user manual are respected. Any use not as prescribed meaning not in accordance with the provisions can cause damage to people or objects. In such cases A.M.I. Italia S.r.l declines all responsibility.

1.3 Guarantee

The device **Geo Saver** has a guarantee of 6 (six)*years.

The non-rechargeable battery Li-SOCl₂ (SAV-C1032) has a guarantee of 4 (four)* years in Stand-by mode (assuming a battery activation test, daily self-tests without turning on the AED). This information refers to new batteries, fully charged at a temperature of 20°C and humidity of 45%.

* For more information consult **Chapter 16** "Guarantee contract for defibrillators Geo Saver Series"

1.4 Exclusion of liability

The rights of liability are excluded in cases of damage to people or objects, if attributable to one of the indicated causes:

- Use of the appliance not in compliance with the provisions.
- Improper use and maintenance of the appliance.
- Use of the device and / or its accessories which show obvious or partial damage.
- Failure to comply with the instructions in the user manual concerning precautions, operation, maintenance and repair of the appliance.
- Use of non-original accessories and/or parts not approved by the manufacturer.
- Arbitrary interventions, repairs or modifications of the device.
- Arbitrary overcoming of performance limits.
- Lack of surveillance of parts subject to wear.

1.5 Indications

The **Geo Saver can** only be used if the patient:

- is unconscious and
- does not breathe and...
- shows no signs of blood circulation

1.6 Counter indications

The **Geo Saver cannot** be used if the patient:

- is in a conscious state or...
- shows normal respiration or...
- shows signs of blood circulation



1.7 Version information

This user manual has a version number. The version number changes every time the manual is updated for changes made to the function of the device or to the device itself. The contents of this user manual are subject to change without notice. The information on the version of this manual is as follows.

Version number: 4.2
Issuing date: 01/09/2020

1.8 Symbols in the manual

In this user manual there are several symbols that indicate the various precautions for use:

SYMBOL	INDICATION	DESCRIPTION
	DANGER	Indicates an immediate risk to the safety of people, which also involves death and damage to the device or parts thereof
	CAUTION	Indicates an unsafe situation or practice involving serious personal injury and damage to the device or parts thereof

1.9 Manufacturer contacts

You can contact our company at the following addresses:

A.M.I. Italia S.r.l.

REGISTERED OFFICE

Via G. Porzio Centro Direzionale Isola G/2
80143 Napoli (NA) Italy

PRODUCTION SITE ITALY
PRODUCTION, LABORATORIES, OFFICES
Via Cupa Reginella, 15/A
80010 Quarto (NA) Italy
Phone: +39 081 806 34 75
Fax: +39 081 876 47 69

PRODUCTION SITE HUNGARY
PRODUCTION, LABORATORIES, OFFICES
A.M.I. International KFT
Közúzó u. 5/A
2000 Szentendre (Hungary)
Phone: +36 26 302.210

Assistance request

email: info@amiitalia.com
Phone: +39 081 806 05 74
Website: www.amiitalia.com

2 Safety instructions

For a correct use of the *Geo Saver* defibrillator, users must be aware of the safety factors listed below.

Please read them carefully.

The *Geo Saver* defibrillator, individually and in connection with its standard and optional (original) accessories, complies with the safety regulations currently in force and is in compliance with the provisions of the directives on medical products.

The appliance and its accessories are to be considered safe in the case of application according to the provisions and if the descriptions and indications listed in this user manual are respected.

The following are the main precautions to be taken for the correct and safe use of the defibrillator, divided for easy consultation between hazard indications, warning indications and disposal instructions.

2.1 Indications of **DANGER**



- Use the *Geo Saver* in accordance with the prescriptions in this user manual. Carefully read these instructions for use and in particular the safety instructions indicated in them.
- In accordance with IEC standards (section 2.4), the use of the *Geo Saver* device or its accessories in the presence of flammable substances (petrol or similar) or in an atmosphere enriched with oxygen or flammable gases / vapours is not allowed.
- Do not recharge the Li-SOCI2 battery (SAV-C1032). Explosion risk!
- Do not allow the batteries to come into contact with an open flame. Do not expose to fire.
- Do not short circuit the battery terminals.
- In case of leakage of liquids or strange smells from the batteries, keep them away from fire to prevent any leaked electrolytes from catching fire.
- Danger of electric shock. The device generates high voltages and dangerous levels of current.
- Do not open the *Geo Saver*, do not remove the panels and do not attempt to repair it. The *Geo Saver* contains no components that users can repair. For repair purposes, the Geo Saver P must be sent to an authorized technical service centre.
- Do not apply the electrodes to the patient's chest if nitro-glycerine patches are present. Remove the patches and only then position the electrodes. Otherwise there is a risk of causing an explosion.
- Do not touch the patient and prevent third parties from coming into contact with the patient during the defibrillation shock phase. Avoid any contact between:
 - parts of the patient's body
 - conductive liquids (such as gel, blood or solution of table salt)
 - metal objects in the surroundings of the patient (such as bed frame or stretching device) that represent indirect ways for the defibrillation current
- Before using the device ensure the patient's safety, if necessary move them carefully and position them in a safe place as per the AHA / ERC 2017 guidelines
- Do not immerse any part of the *Geo Saver*, its parts or accessories in water or other liquids.
- Do not allow liquids to enter the *Geo Saver* its parts or accessories. Avoid spilling liquid on the device and its accessories. Failure to do so may cause damage or cause a risk of fire or electric shock. Do not sterilize the *Geo Saver* or its accessories.



2.2 Indications of **CAUTION**

- Avoid the formation of air bubbles between the skin and defibrillation PADS. The formation of air bubbles during defibrillation can cause severe burns to the patient's epidermis. To avoid the formation of air bubbles, make sure that the electrodes fully adhere to the skin. Do not use electrodes whose gel has dried, check the expiration date before use.
- Do not delay treatment in patients with an implanted pacemaker and perform a defibrillation attempt if the patient has lost consciousness and is not breathing or breathing normally. The *Geo Saver* is equipped with a pacemaker detection system that allows ignoring the signal emitted by the latter; however, with some types of pacemakers, *Geo Saver* may discourage a defibrillation shock.

During the application of the electrodes:

- Do not apply the electrodes directly to an implanted device.
 - Apply the electrodes at least 2.54 cm (1 inch) from any implanted device
- RF (radio frequency) interference, caused by devices such as cellular phones and two-way radios, can cause the *Geo Saver* to malfunction. The *Geo Saver* must be kept at least 2 meters away from these RF devices, as indicated in the standards of EN 61000- 4-3. Keep away from other therapeutic and diagnostic energy sources (eg diathermy, high-frequency surgery, magnetic tomography).
 - Use the *Geo Saver* only if you have achieved a BLS-D or ALS-D training course.
 - Before using the device, make sure that there is no obvious damage.
 - The infrared interface emits optically invisible radiation. The emission diode complies with IEC/EN 60825-1 Class 1 "Eye Save"
 - Do not use paediatric defibrillation PADS (SAV-C0016) on adult patients (older than 8 years and weighing more than 25Kg). Using paediatric defibrillation PADS the *Geo Saver* automatically switches to paediatric mode, reducing the maximum energy available to 50J.
 - Arrange the patient cables so as to reduce the possibility of wrapping or strangling the patient.

- In a domestic environment, keep the defibrillator out of the reach of children and pets.
- Do not apply the defibrillation electrodes directly on an implanted pacemaker to avoid any errors in the interpretation of the device and to avoid damage to the pacemaker through the defibrillation impulse.
- Disconnect high-voltage pulse-sensitive equipment from the patient, ie that is not defibrillator-proof, before delivering the shock.

CAUTION



- Do not allow defibrillation electrodes to touch or come into contact with ECG electrodes, swabs, transdermal patches, etc. Failure to do so may result in creation of electric arcs and burns to the patient during defibrillation, and even current leakage.
- Position the defibrillation PADS as indicated in this user manual and indicated on the package.
- Do not use defibrillation PADS if the gel has been detached from the support or is torn, split or dry.
- If damage has been detected, do not operate the **Geo Saver** under any circumstances.
- Before using the device, remove metal objects from the patient's body (including necklaces or bracelets, etc.)
- Do not use defibrillation PADS other than those supplied by the manufacturer. Otherwise the defibrillator may make false interpretations.
- Do not use defibrillation PADS if they are damaged, even partially.
- Do not touch the patient or PADS during heart rhythm analysis.
- Moving or transporting the patient during the cardiac rhythm analysis performed by the device can lead to an incorrect or not timely diagnosis. During the heart rhythm analysis phase, minimize the movements. If the device is used in an ambulance in motion, stop the vehicle and start again only after having delivered the shock.
- In order to use the **Geo Saver**, you must have completed a training course for basic or advanced cardio-pulmonary resuscitation with the use of a defibrillator (BLS-D or ALS-D course)
- Before applying the defibrillation PADS, if necessary, dry the patient's chest and remove unwanted hair.
- Do not subject **Geo Saver**, its accessories, its parts to falls and / or strong impacts
- Do not use damaged accessories and / or parts, otherwise the device may malfunction.
- Use only original accessories and / or spare parts.
- Avoid excessively aggressive handling of the device of its accessories or parts in order to avoid possible damage. Inspect the entire system periodically.
- Carry out the sanitation operations of the device in compliance with the standards indicated in paragraph 10.3 and always make sure that the device is switched off with the battery removed and PADS disconnected.
- Defibrillation PADS are disposable, to be used only on one patient. Do not reuse defibrillation pads; discard after use and replace with a new pair.
- Defibrillation PADS are not sterile or sterilizable.
- Intense or prolonged administration of cardiopulmonary resuscitation with defibrillation electrodes applied to the patient can damage the electrodes. Replace them if they are damaged due to use or handling.
- Improper maintenance can damage the **Geo Saver** or cause it to malfunction. Follow the instructions in this user manual.
- Use original non-rechargeable Li-SOCI2 (SAV-C1032) batteries from A.M.I. Italia S.r.l. before the indicated expiration date.
- Recharge the rechargeable Li-ion battery (SAV-C1033) at least once every 4 months ensure its perfect function and extend its life.
- The Li-ion rechargeable batteries ACC model (SAV-C1033) must be charged using only the SAV-C1035 battery charger from A.M.I. Italia S.r.l. otherwise the batteries could be damaged
- Remove the batteries from the device only if it has been turned off for at least 5 seconds. Otherwise the device and the battery can be seriously damaged.
- The **Geo Saver**, its parts and accessories are not sterile or sterilizable
- Do not expose the **Geo Saver**, its parts or accessories to direct light or high temperatures
- The Battery Charger (SAV-C1035) must only be used with the Meanwell power supply model GS40A15-P1J (SAV-C1037) supplied by A.M.I. Italia S.r.l. The use of different power supplies could compromise the correct functioning of the battery charger and damage the ACC rechargeable batteries (SAV-C1033)
- In order to safeguard the battery life (SAV-C1032) and guarantee automatic daily tests, after installing it, it is advisable to not remove the battery (SAV-C1032) unless it is to be replaced. The removal of the battery and the subsequent insertion involves a complete test of the AED which considerably consumes its capacity. Furthermore, if the battery is not properly attached it could be damaged.



2.3 Indications of **DISPOSAL**

- The **Geo Saver**, its parts and accessories must not be disposed of with other household waste within the European community. To prevent possible damage to the environment or human health caused by incorrect waste disposal, recycle this product responsibly also to promote sustainable use of resources. To dispose of the used product, use the appropriate waste collection services or return it to the local distributor. In this way it will be possible to recycle safely for the environment

3 Description of the device

3.1 Device Information

Geo Saver is called **AED** or **Automated External Defibrillator** equipped with a mini LCD display.

The device was designed to be used by lay personnel as well as by healthcare personnel who have duly obtained and passed a BLS course according to international guidelines.

Designed to automatically detect and analyse the victim's heart rhythm, it is capable of delivering one or more defibrillation shock if ventricular fibrillation or ventricular tachycardia (monomorphic or polymorphic with beat > 180) is detected. The energy is supplied by an exponential truncated biphasic electrical shock (B.T.E.) able to adapt to the patient's thoracic impedance.

Geo Saver semi-automatic is available in two versions:

Geo Saver 200J (SGS-B0988) – Maximum deliverable energy 200J

Geo Saver 360J (SGS-B0989) – Maximum deliverable energy 360J

It can be used with two types of battery:

- **Non-rechargeable battery Li-SOCl₂** (SAV-C1032), which requires no maintenance, is guaranteed to operate in standby mode for 4 years or carry out a high number of shock
- **ACC Rechargeable battery Li-ion** (SAV-C1033), recommended for those who use the defibrillator intensively

The device allows the data of the rescue to be recorded on an SD Memory Card so that they can be re-displayed on a PC. During the non-use phase, if battery is installed, the device carries out daily self-tests to verify its functional condition, in order to guarantee its prompt use in the moment of need. On the keyboard of the device there is a mini LCD display and a two-colour LED (red / green) through which it is possible to see the outcome of the functional tests and to know the functional status of the device even if switched off (stand-by mode).

The Geo Saver is a defibrillator equipped with a SIM card, GPS and a battery power system. These characteristics make it an extremely versatile device suitable for installation both in public practices and in constantly moving vehicles such as trains or buses.

The presence of a SIM card allows the Geo Saver to transmit and receive data through the mobile telephone network. The GPS makes it possible to track the movements.

The battery power supply of the system dedicated to Geolocation and communication via mobile phone network is autonomous and additional to that of the components dedicated to the basic functionality (defibrillator).

The information that the Geo Saver sends remotely are consultable through the web platform www.AMISAVERCLOUD.com which is accessible from any device with a browser (PC, tablet or smartphone); no dedicated equipment or software is required.

Each user will be able to access the AMISAVERCLOUD platform in a safe way through the use of personalized credentials; then he will be able to quickly visualize the location and status of their Geo Savers.

Compared to a classic defibrillator, the main differences are:

- Remote control
- Telemetry
- Geo-location
- Teleservice
- Management of AEDs and AEDs' users
- Updates through remote connection

Remote control

Thanks to the AMISAVERCLOUD it is possible to control the device quickly and easily.

By accessing the section dedicated to each Geo Saver it is possible to configure each aspect by sending simple commands and knowing the last ones sent.

Telemetry

The Geo Saver connects with the portal at least once a day and sends an AED log containing its parameters and status information which are consultable through the access to the AMISAVERCLOUD:

- Each Geo Saver is represented by a coloured icon which will allow the immediate recognition of its status :
 - Green: ready to use device;
 - Yellow: Warning of an anomaly that does not affect the defibrillator functionality;
 - Red: device failed, service required;
- The consultation of logs allows the user to carry out the check of the device and its accessories.

In case of connection failure (due to lack of network coverage or malfunctions), AMISAVERCLOUD will notify the disservice to the owner via SMS or email.

Geo-location and Anti-theft

For any Geo Saver is possible to:

- locate the last position on a map;
- tracking the movements (auto tracking) and visualize the path on the map through the web portal;
- notify the owner whenever it is moved, by configuring on the AMISAVERCLOUD the “antitheft” functionality. In this case, the web portal will send a message (via SMS or email) as soon as the device detects a movement.

Teleservice

Streaming ECG

The Geo Saver system allows the ECG to be transmitted in real time to AMISAVERCLOUD and, through the web portal, it is possible to monitor the signal.

All the recorded ECGs are saved in the portal and made available for future consultations.

LIVE call

During a rescue it is possible to make a voice call directly from the Geo Saver to the local EMS using a dedicated button.

Through the portal, it is possible to assign up to three telephone numbers to which the call is addressed. The telephone numbers must be set according to the regulations in place in the individual countries where the Geo Saver is installed.

3.2 Service activation procedure

In order for the Geo Saver series defibrillators to fully provide the features for which they were designed, it is necessary to follow the service activation procedure.

1. Check that the SIM card (*) supplied is correctly inserted;
2. Make sure you are in an area covered by a GSM / GPRS signal;
3. Insert the battery for the Geoloc module into the battery compartment;
4. Insert the Geo Saver battery into the battery compartment;
5. Check that the Geo Saver performs the test designed for battery insertion;
6. Turn off the Geo Saver;
7. Wait a minute and check on the AMISAVERCLOUD the arrival of the session log on the page related to the device being activated.

(*)

- *The SIM cards supplied by the manufacturer only work if inserted in Geo Saver type equipment;*
- *Each SIM card supplied by the manufacturer only works if inserted in the Geo Saver to which it is paired;*
- *A SIM card associated with a Geo Saver does not work if inserted in another Geo Saver series device;*

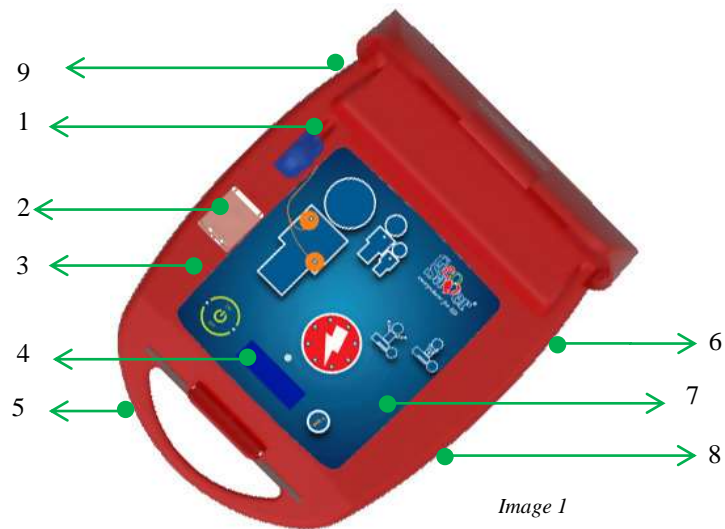
3.3 Classifications

The *Geo Saver* defibrillator is classified as follows:

Code UMDNS	11132
Code GMDN	47910
Code CND	Z12030501
Directory number RDM	1793143 / 1793148
Code CIVAB	DEF01
Class of belonging according to directive 2007/47/CE	II b
Type of protection against electric shock	Internally powered
Type of patient isolation	BF
Degree of protection against penetration of liquids	IPx6
Degree of protection against dust penetration	IP5x
Degree of safety in the presence of a flammable anaesthetic mixture with air, oxygen or nitrous oxide	Not protected
Sterilization or disinfection method suggested by the supplier	See paragraph 11.3
Mode of operation	Continuous operation

4 Description of device details

4.1 General structure of device



Nr.	Description
1	Compartment for PADS connector or ECG cable
2	Live call button
3	Microphone for environmental recordings
4	Status mini display
5	Carrying handle
6	IrDA port (service only, present only in the models with TFT display)
7	Keyboard with buttons and light icons
8	Loudspeaker
9	Live call microphone



Nr.	Description
10	USB port
11	Compartment for SD Memory Card
12	Compartment for SIM card
13	Geo Saver battery compartment
14	Geoloc battery compartment
15	Live call speaker




4.2 Keys, icons and indicators



Image 3

Nr.	Function	Nr.	Function
1	"Patient Type" Indicator Indicates the type of mode in use: Adult if you insert Adult PADS Child if you insert paediatric PADS	6	Product logo Indicates the model of the device
2	"Place PADS" Indicator Icon equipped with LEDs indicating to position the defibrillation PADS.	7	"CPR" Indicator Indicates to start CPR
3	Control LEDs Luminous LED (red / green) allows you to check the functional status of the AED.	8	"Do Not Touch" Indicator Icon equipped with LEDs indicating not to touch the patient during certain operations
4	ON / OFF button Allows you to switch the device on or off.	9	Shock button Equipped with luminous LED allows to deliver a defibrillation shock if indicated
5	Status mini display It allows you to check the functional status of the device	10	"i" button Allows in operating mode to show useful information on the device



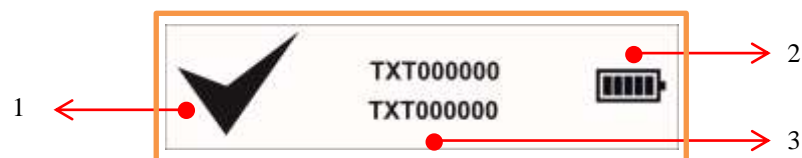



4.3 Status display







Nr.	Description
1	Device status icon
2	Remaining battery level
3	Text commands

4.4 Standard and optional accessories of the device

The *Geo Saver* defibrillator comes with the following standard accessories:

Code	Image	Quantity	Description
SGS-B0988		1 Unit (Version 200J or 360J)	<i>Geo Saver 200J</i>
SGS-B0989			<i>Geo Saver 360J</i>
SAV-C0846		1 Unit	Adult Pads
SAV-C1032		1 Unit	Li-SOCl ₂ AED Disposable battery
SAV-C1073		1 Unit	User guide
SAV-C1038		1 Unit	Li-SOCl ₂ Geoloc disposable battery
SAV-C1071		1 Unit	SIM card

The following are the optional *Geo Saver* accessories that can be purchased separately:

Code	Image	Quantity	Description
SAV-C1033		1 Unit	Li- ion Rechargeable AED battery
SAV-C1035		1 Unit	Charger
SAV-C1037		1 Unit	GS40A15-P1J Power supply
SAV-C1034		1Unit (Contains 3 units)	N.01 Charger
			N.01 GS40A15-P1J Power supply
			N.01 Power supply cable
SAV-C0016		1 Unit	Children Pads
SAV-C0019		1 Unit	CD-ROM Saver View Express
SAV-C0906		1 Unit	SD Card
SAV-C0027		1 Unit	Memory Card reader for PC
SAV-C1072		1 Unit	Carry bag
SAV-C1039		1 Unit	Li-ion rechargeable Geoloc battery
SAV-C1040		1 Unit	Geoloc battery charger

5 Parts and accessories of the Geo Saver

5.1 AED Batteries

The *Geo Saver* defibrillator can work with two different types of batteries:

- (SAV-C1032) Non-rechargeable/Disposable Li-SOCl₂ battery
- (SAV-C1033) Rechargeable Li ion battery

5.1.1 Non rechargeable Li-SOCl₂ battery (SAV-C1032)

The non-rechargeable battery with Li-SOCl₂ technology (SAV-C1032) is supplied fully charged and ready for use. The Li-SOCl₂ non-rechargeable battery has been designed to have a long battery life and no maintenance whatsoever.

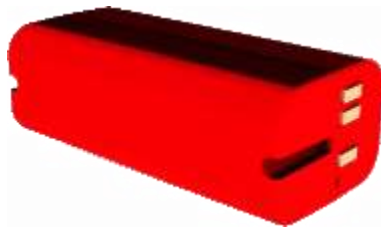


Image 4

The non-rechargeable battery of the *Geo Saver* in Standby mode is guaranteed for 4 (four) years*¹ assuming a battery activation test, daily self-tests without turning on the AED. The Li-SOCl₂ non-rechargeable battery (SAV-C1032) is able to carry out a large number of shocks, which vary according to the version:

<i>Geo Saver Standard 200J</i>	300 complete rescue cycles (shocks at 200J. and CPR)* ¹
<i>Geo Saver Power 360J</i>	200 complete rescue cycles (shocks at 360J. and CPR)* ¹

*¹New and fully charged battery, constant temperature at 20°C and relative humidity without condensation 45%

If the remaining battery level is low, the *Geo Saver* informs the user via audio and visual messages.

The *Geo Saver* will give a **low battery** warning when the level is ≤5% (**WARNING**) and a **very low battery** warning when the level is ≤ 1% (**ALARM**)

WARNING: Remaining capacity level of battery equal or less than 5%.
This notice will only be provided in Operating mode as indicated in paragraph 6.1.
With a 5% battery the *Geo Saver* allows to shock about **14 shocks** or **40 days of stand-by***²

ALARM: Remaining capacity level of Battery at ≤ 1%
This warning will be provided both in Stand-by and in operating mode, as indicated in paragraph 6.1
With a battery at ≤ 1% the *Geo Saver* carries out about **7 shocks** or **20 days of stand-by***²
In this condition the use of the device is not recommended.

*², Constant temperature at 20°C and relative humidity without condensation 45%

!!ATTENTION!!

In order to protect the battery life (SAV-C1032) and guarantee automatic daily tests, after installing it, it is advisable not to remove the battery (SAV-C1032) unless it is replaced. The removal of the battery and the subsequent insertion involves a complete test of the AED, which considerably consumes its capacity. Furthermore, if the battery is not properly attached it could be damaged.

5.1.2 Rechargeable Li-Ion battery (SAV-C1033)

The rechargeable battery with Li-Ion technology (SAV-C1033) of the *Geo Saver* is suitable for those who use the defibrillator intensively. Being rechargeable, it allows operators to reduce management costs and guarantee a greater number of interventions.



Image 5

The ACC rechargeable battery of the *Geo Saver* can be recharged using only the dedicated charger (SAV-C1035) with relative accessories supplied by A.M.I. Italia S.r.l. The battery allows you to carry out a high number of shocks which varies according to the version of the *Geo Saver* in your possession:

<i>Geo Saver Standard 200J</i>	typically 250 continuous shock ^{*1}
<i>Geo Saver Power 360J</i>	typically 160 continuous shock ^{*1}

^{*1}New and fully charged battery, constant temperature at 20°C and relative humidity without condensation 45%

If the remaining battery level is low, the *Geo Saver* informs the user via audio and visual messages. The *Geo Saver* will give a **low battery** warning when the level is $\leq 5\%$ (**WARNING**) and a **very low battery** warning when the level is $\leq 1\%$ (**ALARM**)

WARNING: Remaining capacity level of Battery equal or less than **5%**.
 This notice will only be provided in Operating mode as indicated in paragraph 6.1.
 With a 5% battery the *Geo Saver* allows to shock about **14 shocks** or **40 days of stand-by**^{*2}

ALARM: Remaining capacity level of Battery at $\leq 1\%$
 This warning will be provided both in Stand-by and in operating mode, as indicated in paragraph 6.1
 With a battery at $\leq 1\%$ the *Geo Saver* carries out about **7 shocks/20 days of stand-by**^{*2}
 In this condition the use of the device is not recommended.

^{*2}Constant temperature at 20°C and relative humidity without condensation 45%

It is advisable to replace these batteries every 2 years or after having made a number of recharges greater than **300** (the event that occurs first).

5.1.3 Suggestions for a proper maintenance of battery SAV-C1033

A.M.I Italia recommend that batteries SAV-C1033 left in a “storage stage” to be fully recharged at least every 4 months from the receipt of the goods and to be recharged regularly every 4 months when attached to the device "ready to use" , to avoid completely discharging it and to maintain maximum life expectancy of the battery. The battery pack technology and the modules offered are to ensure a long lasting duration but they require a correct maintenance; failure to follow these requirements will result in an early deterioration of the battery, which will not be covered by warranty.

For warranty replacement consideration, batteries are to be returned to the original supplying distributors/dealer.

5.1.4 Inserting and removing the batteries

To be able to operate the *Geo Saver* the insertion of a battery is required. Below are detailed instructions for correctly installing the batteries (rechargeable or non-rechargeable) in the *Geo Saver*.



Image 6



Image 7



Image 8

- Position the appliance on its side as shown in image (6)
- Remove the cover of the battery compartment as shown in image (6)
- Insert the battery as shown in image (7)
- Push the battery as shown in image (7) positioning it at the bottom of the special compartment
- Close the battery compartment cover as shown in image (8)

Follow the instructions below to **remove** the battery in the device:



Image 9



Image 10

- Position the device as shown in image (9)
- Pull the opposite hooks in order to pull out the battery as shown in image (10)

5.1.5 Recharging station for rechargeable batteries

The charging station (SAV-C1034) allows you to recharge rechargeable batteries with Li-ion technology ACC model (SAV-C1033) of the *Geo Saver*. The charging station consists of the following parts:

- Charger (SAV-C1035) image (11)
- AC/DC power supply/adapter model GS40A15-P1J (SAV-C1037) image (12)
- Power cable with three-pole Italian plug (SAV-C0366) image (13)



Image 11



Image 12



Image 13

5.1.6 Battery charger structure



Image 14

Nr.	Description	Function
1	Charger LED	Indicates the battery charge level or the functional status of battery charger
2	Supply	Input for 12V, 5A power supply connector
3	Battery contacts	Contacts for energy exchange between charger and battery

The Battery Charger (SAV-C1035) must only be used with the AC/DC power supply/adapter supplied by A.M.I. Italia S.r.l. model GS40A15-P1J of Meanwell (SAV-C1037).

The battery charger (SAV-C1035) and the relative power supply unit (SAV-C1037) are not certified under the supervision of the IMQ notified body, therefore they do not fall into the EC certificate no.1104 / MDD. Furthermore, these devices do not have the IMQ mark, therefore they are not indicated in the IMQ certificate no. CA10.00185.

5.1.7 Recharge procedure

- A** Place the charger on a perfectly horizontal shelf and firmly attached to the floor
- B** Connect the power supply (SAV-C1037) to the charger and then to the power outlet
- C** The LED on the charger will flash green, indicating that it is ready to charge
- D** Insert the battery to be charged into the battery charger as shown in image (15)



Image 15

The recharging station allows you to recharge exclusively original ACC rechargeable Li ion batteries (SAV-C1033) of A.M.I. Italia S.r.l. The charging time of around 2.5 hours may increase in the case of batteries that have undergone recharging cycles higher than the one indicated. The battery charger is equipped with a control LED that indicates both its functional status and the battery charge level, if inserted. The following is a table that allows identification of the control LED coding:

INDICATOR	RED		GREEN	
FIXED	Battery not working		Battery charge completed	
FLASHING	Battery inserted	Faulty battery charger	Battery inserted	Battery charging
	Battery not inserted		Battery not inserted	Battery charger waiting for battery insertion

When recharging, the battery charger control LED will flash green with a different frequency depending on the level of recharge, until the charge is fully indicated by the control LED with FIXED green light.

Charge level	0%	25%	50%	75%	100%
Number of consecutive flashes	1	2	3	4	Fixed

5.2 Geoloc Batteries

The *Geoloc* Geo Locator can work with two different types of batteries:

- (SAV-C1038) Non-rechargeable Li-SOCl₂ battery
- (SAV-C1039) ACC Rechargeable Li ion battery

5.2.1 Non-rechargeable Li-SOCl₂ battery (SAV-C1038)

The non-rechargeable battery with Li-SOCl₂ technology (SAV-C1038) is supplied fully charged and ready for use. The Li-SOCl₂ non-rechargeable battery has been designed to have a long battery life and no maintenance whatsoever.



Image 16

The *Geoloc* non-rechargeable battery in Standby mode is guaranteed for 4 (four) years*¹ assuming a machine activation test, daily self-test without any connection to the AMISAVERCLOUD for transfers other than the session log and a good coverage area for the GPS and GPRS signal.

If the remaining battery level is low, the AMISAVERCLOUD highlights the status by colouring yellow or red the graphic representation (pie chart) of the remaining capacity and eventually informs the user by sending e-mail or sms if sending was enabled from the specific configuration panel (par. 10.4).

5.2.2 Rechargeable Li ion battery (SAV-C1039)

The rechargeable battery with Li-ion technology (SAV-C1039) of the *Geoloc* is recommended for those who intensively use the unique features of the Geo Saver- Geoloc- AMISAVERCLOUD (Streaming, Speakerphone calls, tracking). Being rechargeable, it allows operators to reduce management costs and guarantees a greater number of connections.



Image 17

The ACC rechargeable battery of the *Geo Saver* can be recharged using only the dedicated charger (SAV-C1040) with relative accessories supplied by A.M.I. Italia S.r.l.

If the remaining battery level is low, the AMISAVERCLOUD highlights the status by colouring yellow or red the graphic representation (pie chart) of the remaining capacity and eventually informs the user by sending e-mail or sms if sending was enabled from the specific configuration panel (par. 10.4).

!!ATTENTION!!

The rechargeable batteries require maintenance, which consists of a complete recharge every 4 months

5.2.3 Inserting and removing the batteries

In order to function, the *Geoloc* module requires the insertion of a battery. Below are the detailed instructions for correctly installing the batteries (rechargeable or non-rechargeable) for the *Geoloc* module inside the *Geo Saver*.



Image 18



Image 19

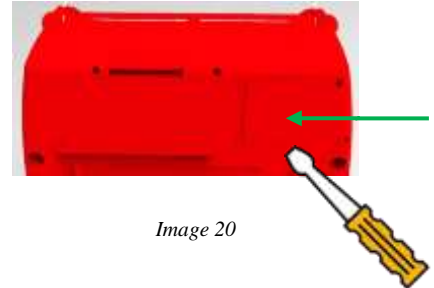


Image 20

- Position the appliance on its side as shown in image(18)
- Remove the cover from the battery compartment as shown in image (18)
- Hook the Battery connector to the Geoloc as shown in image (19)
- Insert the battery
- Close the battery compartment cover as shown in image (20)

Follow the instructions below to remove the battery from the device:



Image 21



Image 22

- Open the cover of the Geoloc battery compartment as shown in image (21)
- Remove the battery
- Disconnect the battery connector as shown in image (22)

5.2.4 Recharge procedure

- A Remove the rechargeable battery (SAV-C1039) from the Geo Saver
- B Connect the battery charger (SAV-C1040) to the power outlet
- C Hook the connector of the battery charger / power supply into the battery as shown in image (23)

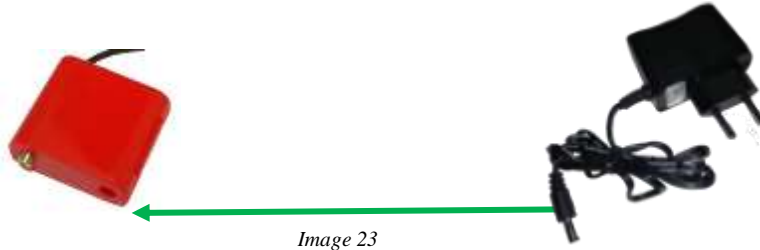


Image 23

The recharging station allows recharging exclusively original ACC rechargeable Li ion batteries (SAV-C1039) of A.M.I. Italia S.r.l. The charging time of around 2.5 hours may increase in the case of batteries that have undergone recharging cycles higher than the one indicated. The battery charger is equipped with a control LED that indicates both its functional status and the battery charge level, if inserted. The following is a diagram that allows the identification of the control LED coding:

INDICATOR	RED		GREEN	
FIXED	Battery inserted	Battery charging	Battery inserted	Battery end of charge
	Battery not inserted		Battery charger waiting for battery insertion	

5.3 Defibrillation PADS

The *Geo Saver* allows the use of two different defibrillation PADS depending on the patient to be treated:

- Defibrillation PADS for **Adults** model SAV-C0846
- Defibrillation PADS for **Children** model SAV-C0016

5.3.1 Defibrillation PADS for Adults SAV-C0846

The SAV-C0846 defibrillation PADS are pre-gelled disposable types.

They must be used on adult patients (**age > 8 years or weight > 25Kg**). Defibrillation PADS are supplied in a single sealed package with the expiration date (typically 30 months). On the expiry date the PADS must be replaced even if not used.

The **SAV-C0846** Pads are characterized by the presence of the cable and the PAD connector **outside the sealed package**. This solution has been adopted in order to maximally speed up the positioning of the Pads avoiding the need to insert the connector during the phases of the rescue.



The **SAV-C0846** PADS are polarized type, that is, the positioning of the electrodes **must not be reversed**. For more information on positioning the defibrillation PADS, refer to the relevant paragraph.

5.3.2 PADS for Children SAV-C0016

The SAV-C0016 defibrillation PADS are pre-gelled disposable types.

They must only be used on children patients (**age < 8 years or weight < 25Kg**). The defibrillation PADS are supplied in a single sealed package with the expiration date (typically 30 months). On the expiry date the PADS must be replaced even if not used.

The cable, connector and PAD adapter **are inside the sealed package**.



The PADS SAV-C0016 allow to deliver shocks on paediatric patients with a maximum energy level of 50J as prescribed by the international guidelines ERC/AHA 2017. The PADS SAV-C0016 are polarized type, that is, the positioning of the electrodes **must not be reversed**.

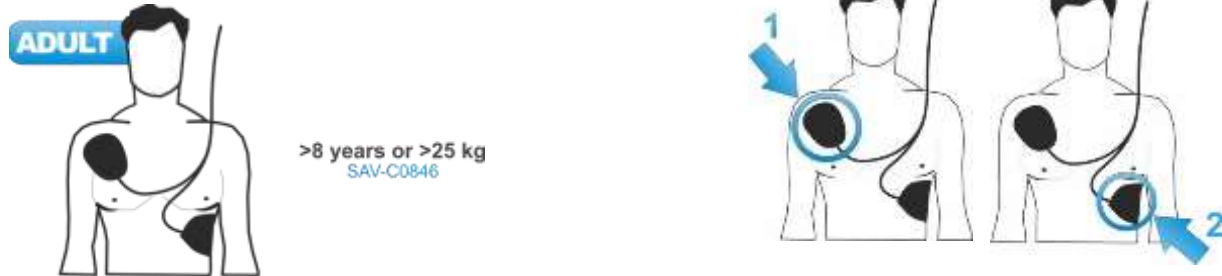
For more information on positioning the defibrillation PADS, refer to the relevant paragraph.

5.3.3 Positioning of defibrillation PADS


The correct placement of the PADS is essential for an efficient analysis of the patient's heart rhythm and for the consequent delivery of the shock (if necessary).

Always refer to the instructions given both on the packaging of the PADS and directly on each individual pad.

The PADS of the *Geo Saver* are polarized type, do not reverse the positioning of each single pad.



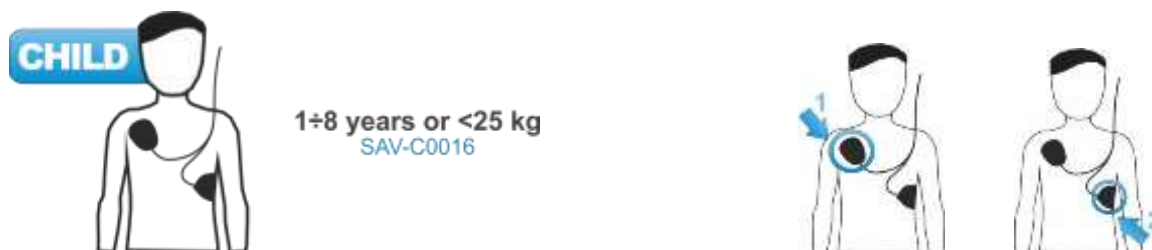- 1 Position **Pad 1** immediately below the patient's right collarbone
- 2 Position **Pad 2** on the ribs on the left side of the patient under the left side of the chest



- 1 Position **Pad 1** immediately below the patient's right collarbone
- 2 Position **Pad 2** on the ribs on the left side of the patient under the left side of the chest

5.4 Memory Card and USB

The *Geo Saver* allows the recording of data on the **internal memory** as well as on the **external memory card**. Supported memory cards are SD/SDHC cards with capacities up to 8GB



Image 24

To install a **Memory Card** in the *Geo Saver* follow this procedure:

- A. The memory card must be inserted before attaching the battery
- B. Place the device on a firm, stable horizontal shelf as shown in the image



Image 25

- C. Insert the Memory Card with the contacts facing upwards as shown in the image and push it until it is completely inserted

The data recorded directly on the internal memory of the *Geo Saver* can be downloaded via the **USB port** on the back of the device (image 25).

The USB cable to be used must be a **mini USB 2.0** (USB / Mini USB connection)



Image 26

To connect the mini USB cables to the *Geo Saver* follow this procedure:

- A. Detach the battery and insert the Mini USB terminal of the cable in the appropriate compartment on the *Geo Saver*
- B. Connect the USB terminal of the cable to a Personal Computer
- C. Use the PC Saver View Express software

!!ATTENTION!!

The USB is a service port used for device configuration purposes (for the exclusive use of personnel authorized by AMI) or for downloading recorded data stored internally.

The functions related to the USB port are enabled and accessible only when the device is turned off.

WARNING: When the device is switched on, and especially when a patient is connected to it, it is recommended:

- do not use the USB port
- do not touch the USB port
- remove the USB cable, if inserted in the USB port, before starting the device

6 Auto test

The Geo Saver has been designed to be a totally safe device, always ready for use and able to automatically and constantly verify correct operation, minimizing maintenance operations.

The Geo Saver performs different types of self-tests:

- **Activation:** Every time a battery is inserted in the device
- **Automatically:** During Stand-by mode with daily/monthly/half-yearly intervals
- **Switching on:** When the device is switched on





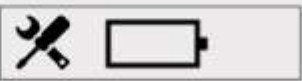








The outcome of the control test can be viewed via a two-colour LED (green/red) and the LCD mini display. The mini-display and the control LED let you know at any time, even when the device is switched off (stand-by mode), the functional status of the device and its battery.

6.1 Display and control LED

Both the display and the control LED are positioned on the front of the **Geo Saver** keyboard.

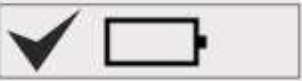Based on the different colour of the control LED and the information shown on the display, the operator can independently determine the functional status of the defibrillator and its battery.

The following table shows the flashing code of the control LED and the relative screens of the control display.

		Control LED green / red blinking	
IN STAND-BY (with battery)	Device ready for use		
	Warning for a low battery level, replace it	 	
	Faulty device, needs service		
WORKING	Device working	OFF	
	Warning: battery is getting low (5% left)	OFF 	 battery is getting low
	Caution! low battery, replace it	 	 low battery, replace it

6.2 ACTIVATION Test


The **Geo Saver** performs functional tests only if the battery is installed.

Each time a battery is inserted, the device will perform a diagnostic ACTIVATION test.

During this test the device performs a complete control (firmware/hardware), which involves a consumption of the battery equal to a shock, therefore it is advisable, once performed, not to remove the battery from the device.

The ACTIVATION test requires a manual intervention by the operator, who must perform the following steps:

❖ Insert the battery into the device

If the battery is correctly inserted, the **Geo Saver** will automatically turn on emitting an acoustic signal and the power button  will light up green while the control LED will turn off.

The following screens will appear on the control display:



The device will emit a voice command (audio) suggesting the operation to complete.

Vocal Command:



“Device test”

“Press the flashing red button”

The shock button will light up with flashing light.



The operator will have a maximum time of 60 seconds* to press the shock button.

❖ Press the shock button to start the ACTIVATION test

If the shock button is pressed correctly it will stop flashing and the device will start the activation test:



If no error is detected the following screen will be shown on the display:



❖ Turn off the device

If it is not to be used immediately turn off the **Geo Saver** and leave the battery inserted to ensure that periodic self-diagnostic tests are performed (see Section 6.3)

If the shock button is not pressed within the time limit indicated by the countdown, the **Geo Saver detects an error.*

Switch the device back on and press the shock button within the time limit indicated by the countdown.

If, on the other hand, the shock button has been pressed but the shock button continues to flash it means that the shock button does not work properly. Turn off the device and perform the operation again; if the problem persists, contact the authorized technical assistance centre.

6.3 AUTOMATIC Test

The *Geo Saver* was designed to always be ready in the moment of real need.

The device is equipped with a mode, called stand-by, in which when the device is switched off (with battery installed) it performs diagnostic tests **automatically** on a daily basis.

The automatic self-test does not require any manual operation by the operator and takes a few seconds.

The *Geo Saver* will inform the operator of the start of the automatic self-test through the mini Control Display:



The automatic self-test involves a reasonable consumption of the battery.

Since a daily test with complete analysis would lead to excessive battery consumption, three levels of automatic tests have been set: **basic** (daily), **in-depth** (monthly), **complete** (half-yearly).

The result of the automatic self-test can be verified using the LED and the mini-control display located on the device keyboard.


Consult the Led table and the mini-control display shown in paragraph 6.1

6.4 POWER ON Test

The *Geo Saver* performs self-diagnostic tests each time it is turned on.

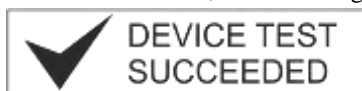
This test is performed in order to verify the proper function of the device before use.

The test is conducted automatically and lasts a few seconds.

After pressing the power button  the *Geo Saver* will emit an acoustic signal to confirm power-up, the control LED will be off and the following screen will be displayed on the colour display:



If no error is found, the following screen will appear on the display:



From this moment the device will be ready to be used and will provide the operator with the first instructions to start the intervention.

If not to be used immediately, switch off the *Geo Saver* and leave the battery inserted to ensure periodic self-testing (see Section 6.3).

7 Information button

The Geo Saver is equipped with an “i” button, through which the operator can view various useful information about the device in use on the LCD display.

This button can only be used when the device is switched on and is automatically disabled in the event of an on-going emergency.

The information shown on the display is divided into three different pages, which can be consulted by pressing the “i” button n times (n meaning the number of pages).

Below is a detailed description of the procedure for using this button and the information displayed:

➤ Switch on the device

The Geo Saver will carry out the automatic power on test after which it will be ready for use.



➤ Press the “i” button

1. After pressing the button the **first time** the following screen will appear with the related information:



MODEL: GEOSAVER 200J
S/N: XXXXXXXXXXXXXXXXXXXX
POWER: BATTERY

Nr.	Description
1	Device model
	Device Serial Number
2	Type of supply

2. Pressing the button for the **second time** will display the following screen with the related information:



PROTOCOL: 150-200-200J
SHOCKS: 6
DATE: 01/03/2019

Nr.	Description
3	Pre-set shock protocol
4	Number of shocks made
5	Current date

3. Pressing the button for the **third time** will display the following screen with the related information:



LANGUAGE -> ITALIAN
Italian
English

Nr.	Description
6	Change language

To change the language, press the "i" button for about 3 seconds and release it.

The user will see the following screen:

LANGUAGE -> ITALIAN
-> Italian
English

Select the desired language by pressing the "i" button. Once the desired language has been selected, press the "i" button for about 3 seconds to confirm the selection.

The selected language remains in memory and at the next restart of the device it will be the default language.

8 Defibrillation

If you need to help a person with sudden cardiac arrest, remember to implement the sequence of actions recommended by ERC and AHA 2017.

The ERC has established a rescue protocol to be respected during the resuscitation of a person suffering from sudden cardiac arrest. This protocol has been called the "chain of life"



Image 27

- 1 Make sure the person needs help (*unconscious and not breathing and does not shows signs of blood circulation*)and call 118 immediately
- 2 While waiting for a defibrillator to be available, begin CPR manoeuvres immediately
- 3 Use the **Geo Saver** defibrillator to restore normal heart rhythm
- 4 Continue until resuscitation of medical competence

8.1 Switching on the Geo Saver

The **Geo Saver** defibrillator will automatically start the semi-automatic defibrillation mode every time it is turned on (default setting). Below are the procedures to follow to use the device in this mode.

Press the power button on the device

The **Geo Saver** will emit an acoustic signal to confirm the power on, the ON / OFF button will be lit fixed green. The following screens will appear on the LCD display:



If the test is successful, the device will suggest the first operation to be performed by the operator using voice (audio) and visual commands (luminous icons and display text):

Voice commands	Display	Luminous Icons Keyboard
Make the emergency call	call now the emergency service	Command Position Defibrillator PADS
Remain calm and follow the voice instructions. If the patient is unconscious and does not breathe, remove clothing in order to apply the electrodes on the bare chest	Remove clothing to expose bare chest	
Place the electrodes firmly on the chest as shown in the figure	Place Electrodes firmly to bare chest	

8.2 Preparation of the patient

In order to correctly position the defibrillation PADS on the chest, it is necessary to carry out the following preliminary operations:

Remove clothing from the patient's chest

If the patient's chest has thick hair it might be necessary to shave the Pads positioning points.


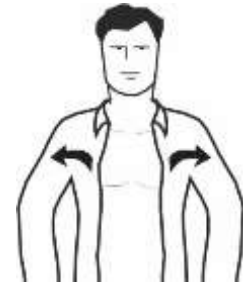



Image 28

8.3 Positioning the PADS

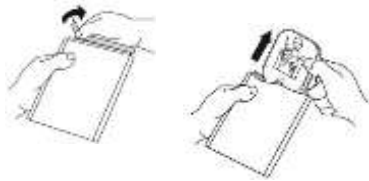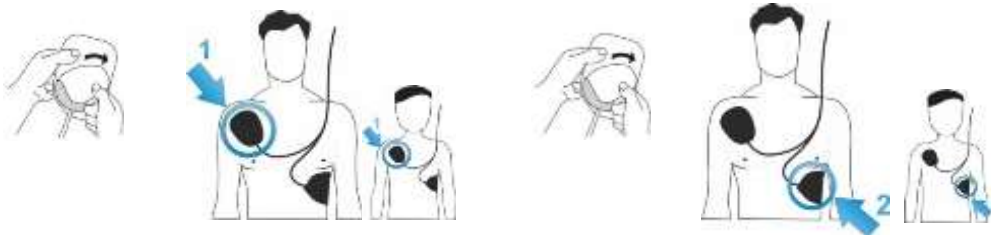A Remove the defibrillation PADS from their packaging.




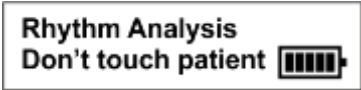
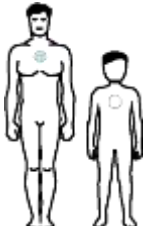
B Remove the protective film from each individual PAD and place it on the patient's chest



The defibrillation PADS are of the polarized type, this requires that they be positioned at the points indicated on each single PAD. If the patient is a child, use the SAV-C0016 paediatric defibrillation PADS, see the relative paragraph for more information.

8.4 Cardiac rhythm analysis

The *Geo Saver* automatically analyses the patient's heart rhythm. This phase is highlighted by the following commands:

Voice commands	Display	Luminous Icons Keyboard
Do not touch the patient		 <p>Icon do not touch the patient on with fixed light</p>
Cardiac rhythm analysis in progress		 <p>Icon type of pads used (<i>adult / child</i>)</p>

During cardiac rhythm analysis the patient's body must not be touched and must not be subjected to vibrations or movements. The *Geo Saver* analysis software was designed to recommend defibrillation shock treatment only if the patient is suffering from the following arrhythmias:

VF *Ventricular Fibrillation*

Peak to peak Amplitude min. 200 μ Volts



Some rhythms with very low amplitude or low frequency VF may not be interpreted as shockable.

VT *Ventricular tachycardia* (including ventricular flutter and polymorphic ventricular tachycardia)

Rhythm frequency min. 180 bpm and peak-to-peak amplitude min. 200 μ Volts



Some rhythms with very low amplitude or low frequency VT may not be interpreted as shockable.







The *Geo Saver* detects both noise artefacts in the ECG, caused, for example, by patient movement, defibrillation electrode adjustment, and electronic disturbances from external sources. In these cases the analysis is delayed or interrupted.

The *Geo Saver* is able to detect and filter impulses coming from an implanted pacemaker.

8.5 Shockable Rhythm

If a Ventricular Fibrillation or Tachycardia is detected, the Geo Saver will inform the operator using the following commands:


Voice commands	Display	Icons/ Luminous buttons
Shock recommended	<div style="border: 1px solid black; padding: 5px; display: inline-block;"> device is charging stay clear patient  </div>	 <p>Do not touch the patient icon on fixed light</p>
Keep distance, charging		
Press shock button	<div style="border: 1px solid black; padding: 5px; display: inline-block;"> Press shock Button  </div>	 <p>Shock button flashing in red</p>

To deliver the shock press the shock button within 15 seconds.


Before pressing the shock button move away and make sure no one is touching the patient.

If, the shock button is not pressed within 15 seconds from the shock warning, the *Geo Saver* will automatically disarm.

At this phase the shock button will no longer flash and the device will inform the user through the following commands:

Voice commands	Display
Shock cancelled	<div style="border: 1px solid black; padding: 5px; display: inline-block;"> SHOCK CANCELLED  </div>
Shock button was not pressed	

If, on the other hand, the shock button was pressed the *Geo Saver* will deliver the defibrillation shock. In this phase the shock button will no longer flash and the device will inform the user through the following commands:

Voice commands	Display
Shock delivered	<div style="border: 1px solid black; padding: 5px; display: inline-block;"> SHOCK DELIVERED  </div>

The *Geo Saver* delivers the shock using the BTE waveform with auto compensation of the patient's thoracic impedance. The *Geo Saver* download protocol is incremental, ie the energy delivered to the patient varies incrementally based on the number of shocks performed:

Geo Saver 200J: The first shock is delivered to energy of **150J** the following at **200J**

Geo Saver 360J: The first shock is delivered to energy of **200J** the second at **250J** the following at **360J**


The detected impedance value must be between 20 and 200 Ohm; if a value outside this range is detected, it is required to position the PADS.

The shock protocol is pre-set, cannot be modified by the user and is reset at each power up.

8.6 Change of rhythm

The *Geo Saver* analyses the patient's heart rhythm continuously, during the resuscitation phase.

If the device after recommending the shock detects a change in the heart rate of the patient who no longer needs a defibrillation, it will carry out the automatic disarmament. This phase is highlighted by the following commands:

Voice commands	Display
Shock cancelled	
Rhythm changed	

8.7 Non-shockable rhythm

If the *Geo Saver* during cardiac rhythm analysis does not detect a VF or a VT, it will inform the user through the following commands:

Voice commands
Shock not recommended

All rhythms other than VF and VT will be considered as non-shockable. For more information, paragraph 10.6.

8.8 CPR

The *Geo Saver* defibrillator will guide the operator to CPR (Cardio Pulmonary Resuscitation) in one of the following cases:











- A shockable rhythm was detected and a defibrillation shock was delivered
- A non-shockable rhythm was found
- A shockable rhythm was found, but the shock button was not pressed
- A shockable rhythm was found but the patient's rhythm changed

The *Geo Saver* will provide instant by instant instructions for performing CPR, instructing the operator on how to perform chest compressions and insufflations. According to the 2017 AHA /ERC guidelines, the duration of cardio-pulmonary resuscitation is about 2 minutes.

If resuscitation is carried out by a single operator, the compressions/insufflations ratio must be 30/2 for 2 minutes (5 cycles) for both adult and child patients.

During the phase of chest compressions the *Geo Saver* mediates a metronome which will set the rhythm to maintain in order to perform the compressions at the right time. Once the compressions have ended, you will need to perform the two breaths. These instructions are repeated throughout the CPR phase, or about 2 minutes.

The following table shows the main operations to be performed and the relative visual-text-vocal commands provided by *Geo Saver*

No.	Type of command (<i>Geo Saver</i>)	Instruction <i>Geo Saver</i>	Operations to be performed
1	Vocal	“Begin Cardio-Pulmonary Resuscitation”	<p>A. Make sure the patient is lying on a flat surface B. Kneel beside the victim C. Place the heel of one hand in the center of the victim's chest D. Place the heel of the other hand over the first one E. Interlace the fingers of both hands and make sure that the pressure is not applied to the ribs. Do not exert any pressure on the upper portion of the abdomen or the lower portion of the sternum</p> 
	Visual DISPLAY		
	Visual LUMINOUS ICON		
2	Vocal	“Quickly compress the patient's chest”	<p>F. Stand vertically on the victim's chest and, with arms extended, compress the sternum. Keeping the arms stretched, the external cardiac massage is exercised using the weight of the trunk; the oscillation movement must be from pivoting on the coxo-femoral joint G. After each compression release all pressure from the chest without losing contact between one's hands and the sternum; repeat the maneuver with a frequency of 100 / min (a little less than 2 compressions per second) H. The compression and release phase must be approximately equal in duration</p> 
	Visual DISPLAY		
	Visual LUMINOUS ICON		
	Acoustic Signal (BEEP)	The <i>Geo Saver</i> signals with a BEEP every compression to be performed.	
3	Vocal	“Perform two breaths” “Exhale” “Exhale”	<p>Immediately open the air passage using the head and chin towards the back maneuver</p>  <p>Perform two insufflations The rescuer inhales normally and keeping the chin lifted with two fingers, make the lips adhere around the mouth of the injured person. The contralateral hand closes the nostrils to avoid air release and keeps the head in hyperextension. Blow out the air by performing a normal expiration lasting about 1 second.</p> 
	Visual DISPLAY		
	Visual LUMINOUS ICON		
4	The <i>Geo Saver</i> will repeat STEP 1 to 3 for about 2 minutes		Follow the voice and text instructions of the <i>Geo Saver</i> until the device stops the CPR phase (about 2 minutes)

9 Recording, viewing and archiving the data

The *Geo Saver* defibrillator is able to record and store both **the SERVICE data of the device** and **the complete data of the rescue operations** carried out. Data recording and archiving is done automatically (cannot be not deactivated by the user) both on **the internal memory** of the device and on **the memory card** when installed.

9.1 Data recording

The **internal memory** of the *Geo Saver* allows the storage of up to 6 hours of environmental recordings (audio), ECG tracing, patient data (FC and Ω) and all rescue events. The stored data can be viewed on a PC using the PC Saver View Express software (SAV-C0019).

Two types of files are stored on the **external memory SD Card**:

- **AED1LOG.txt** This file stores all the automatic self-tests performed by the device with its outcomes and all the **SERVICE** information. This type of file can be viewed on a PC using a simple reading program.
- **AEDFILE.aed** This file stores the rescue data such as: environmental recordings (audio), ECG tracing, patient data (FC and Ω) and all rescue events. This type of file can be redisplayed on a PC using the PC Saver View Express software.

The number and duration of recordings depend on the capacity of the memory card, below is an example:

Type	Capacity	Stored Data	
SD Card	512 MB	Sounds, Events, Parameters, ECG. Service (AED1LOG + AEDFILE)	1.500 minutes (25 hours)
	1 GB		3.000 minutes (50 hours)
	2 GB		6.000 minutes (100 hours)
SDHC Card	4 GB		12.000 minutes (200 hours)

9.2 Archiving data on PC

The rescue data recorded by the *Geo Saver* defibrillator can be stored, analysed and printed from a Personal Computer using the management software Saver View Express.



Saver View Express®

Image 29

For more details on the PC Saver View Express software, consult the relevant user manual.

10 Web Management

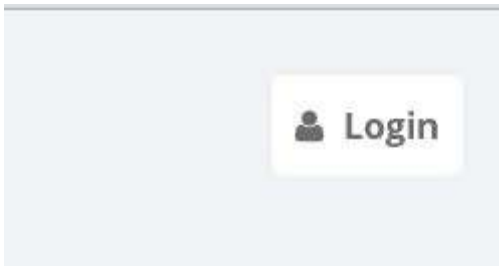
AMISAVERCLOUD is the platform through which a user can check the status of its Geo Saver .

To access the functions of AMISAVERCLOUD there is a mechanism of authentication with password.


When a Geo Saver is purchased, a new user is generated in AMISAVERCLOUD by referring to an e-mail provided by the customer. When the new account is created, the system automatically sends the login credentials with an invitation to change the password at the time of the first login to secure the data relating to the Geo Saver .

To access your account, go to the page:

www.amisavercloud.com



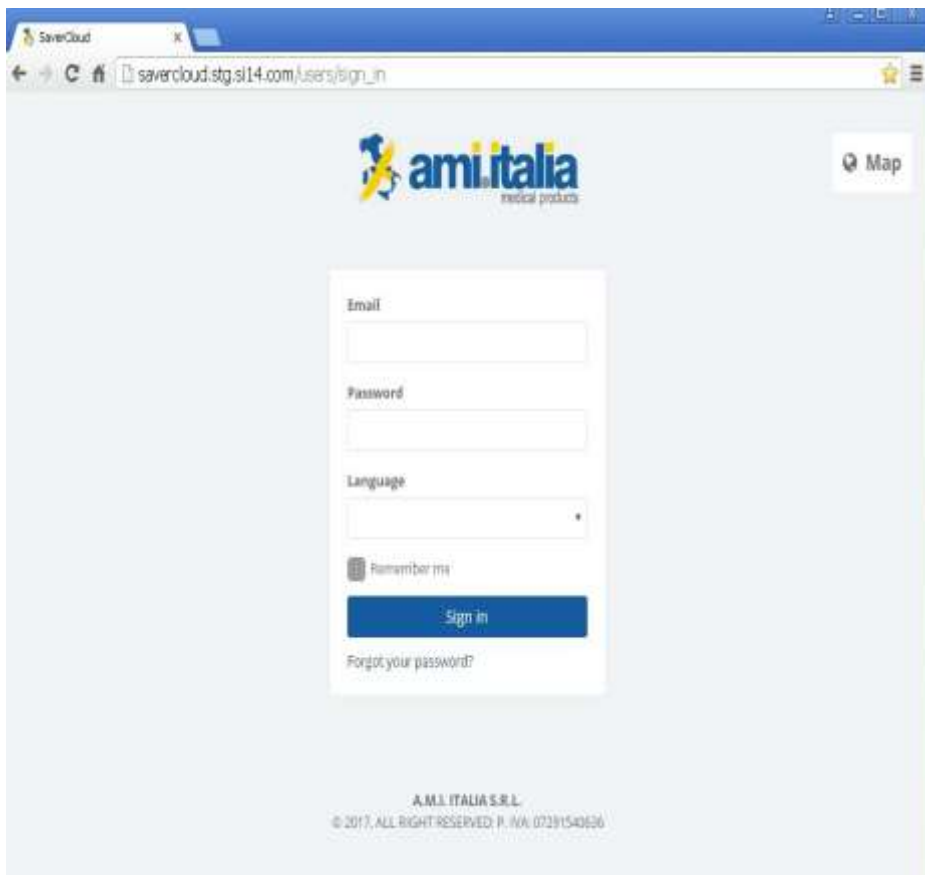Click on the button “Login” and enter your credentials to enter your account.



10.1 Main features

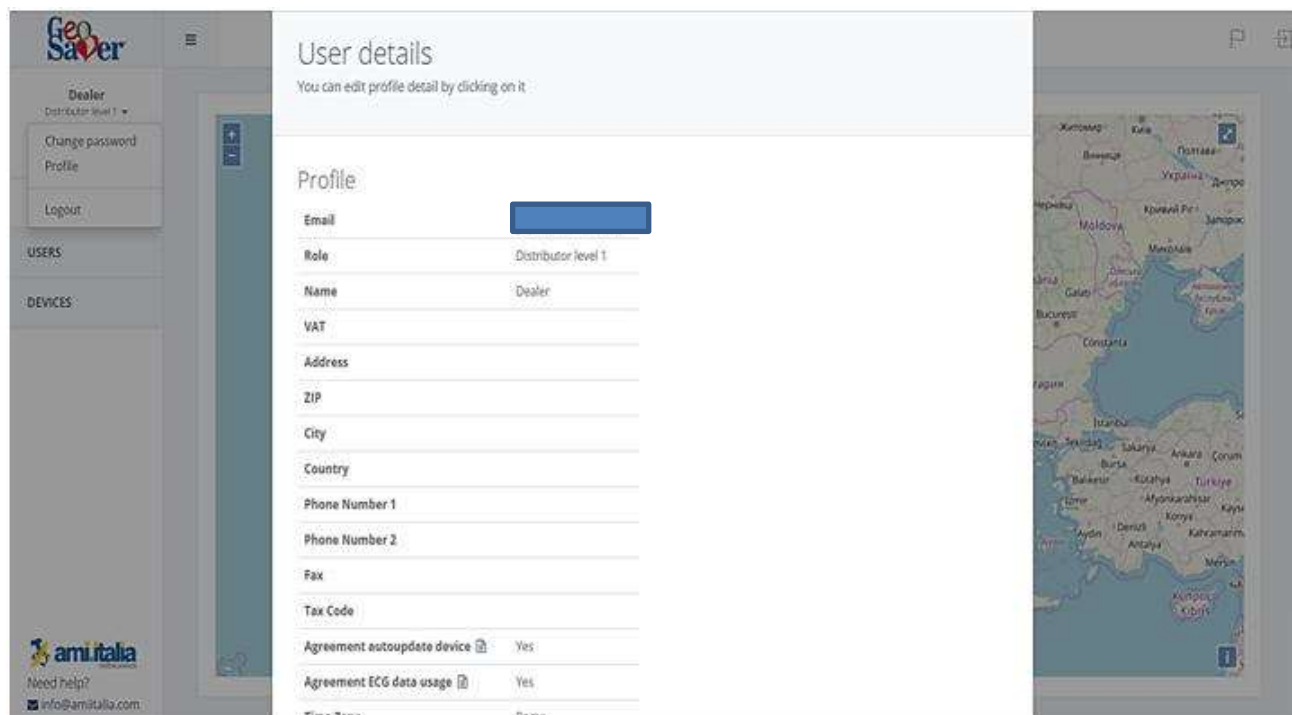
10.1.1 Changing the language

AMISAVERCLOUD is a multilingual cloud service. You can select one of the languages available through the specific icon on the top right




10.1.2 User profile and password

The current user has the possibility to view his profile and modify the password via the links located on the menu.



10.1.3 Device management

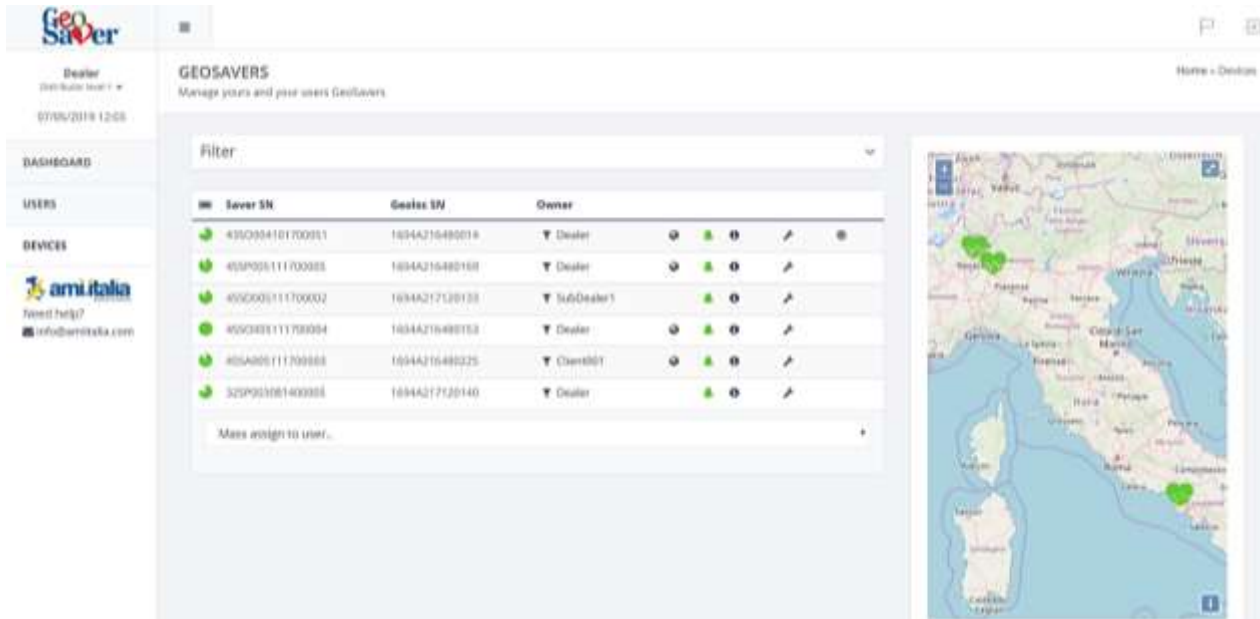Each user can be associated with one or more Geo Savers. These are listed under the "DEVICES" menu item and displayed on a map that briefly summarizes their status.



It is possible to view the details of each device by following the link in the table (or inside the map) for each of them.

When viewing a Geo Saver detail, the menu is expanded with the following items:

- Geo Saver: summary page details
- Saver Log: list of logs received from the Saver
- Geoloc Log: list of logs received from the Geoloc
- ECG Streaming List (if data is present): list of ECG streaming files and possibility to download or reproduce each of them.
- GPS tracks (if data is present): list of GPS tracks recorded by Geo Saver and display of each of these on a map
- Maintenance: management of accessories duration (Geoloc battery, pads)
- Alarms configuration: configuration of the recipients to whom the alarm messages are to be sent (via SMS or e-mail).



10.1.4 ECG streaming


To activate the ECG transmission, switch on the Geo Saver and then enable ECG transmission.

While the ECG is streaming from the Geo Saver, a notification is displayed on the cloud and the "STREAMING" section becomes accessible as a "GEO SAVER" sub-level of the menu.

The device in the ECG check phase is highlighted in the Geo Saver table and a link to the section appears on the same page. In this section it is possible to view the streaming in real time.

At the end of the ECG control phase the transmitted data are saved and becomes available for review in the " ECG STREAM HISTORY " section.

Filename	Acquired on	Duration
2019-04-23_174331_43SO004101700051.stream	April 23, 2019 17:45	00:03:35.396
2019-03-26_180542_43SO004101700051.stream	March 26, 2019 18:10	00:01:02.963
2019-02-21_145848_43SO004101700051.stream	February 21, 2019 14:59	00:01:03.398
2018-11-09_122751_43SO004101700051.stream	November 08, 2018 12:29	00:00:39.878
2018-11-08_164536_43SO004101700051.stream	November 08, 2018 16:46	00:00:43.853

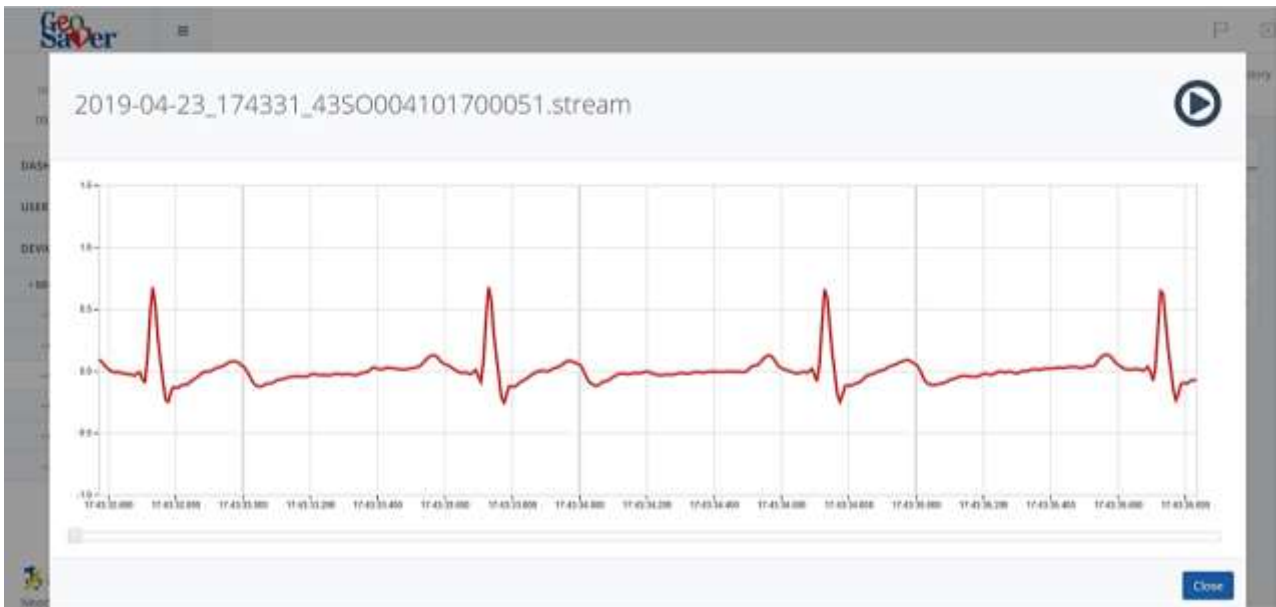


10.1.5 Auto tracking

Among the command that can be sent to the Geo Saver through the cloud there's enable/disable Auto tracking.

This functionality can be found in the GEOSAVER page –GEOLOC section.


When this function is active the Geo Saver, after having detected a movement, will start sending its positions to the cloud. After few minutes from the last move, the file with the path becomes available in the "GPS TRACKS" section accessible from the Geo Saver page as a menu sub-level.

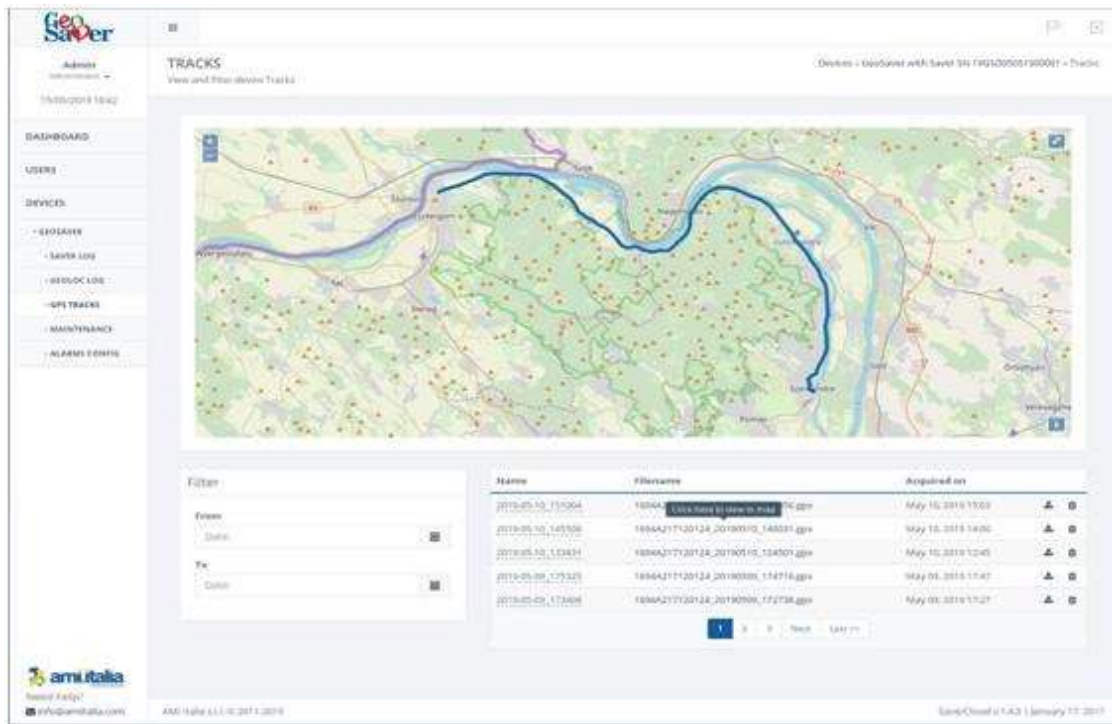Clicking on one of the files in the table shows the path on the map.



10.1.6 Battery replacement

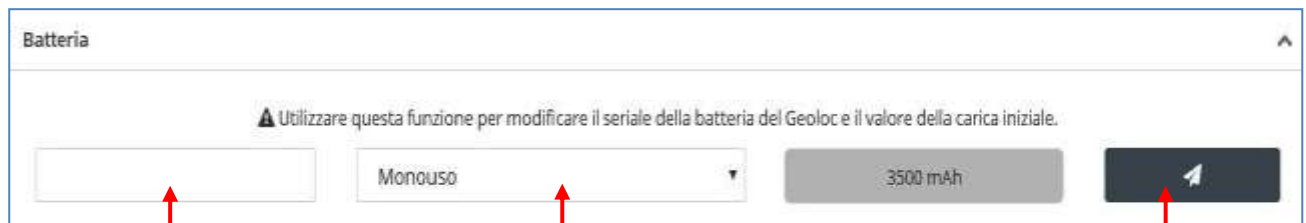Follow the instructions in paragraph 5.2.3 to replace the Geoloc battery. Once completed, go to the AMISAVERCLOUD and, from the MAINTENANCE page of the device in possession, follow the instructions below.



1

2

3

Enter the Serial Number of the inserted battery here

Select the type of battery inserted here

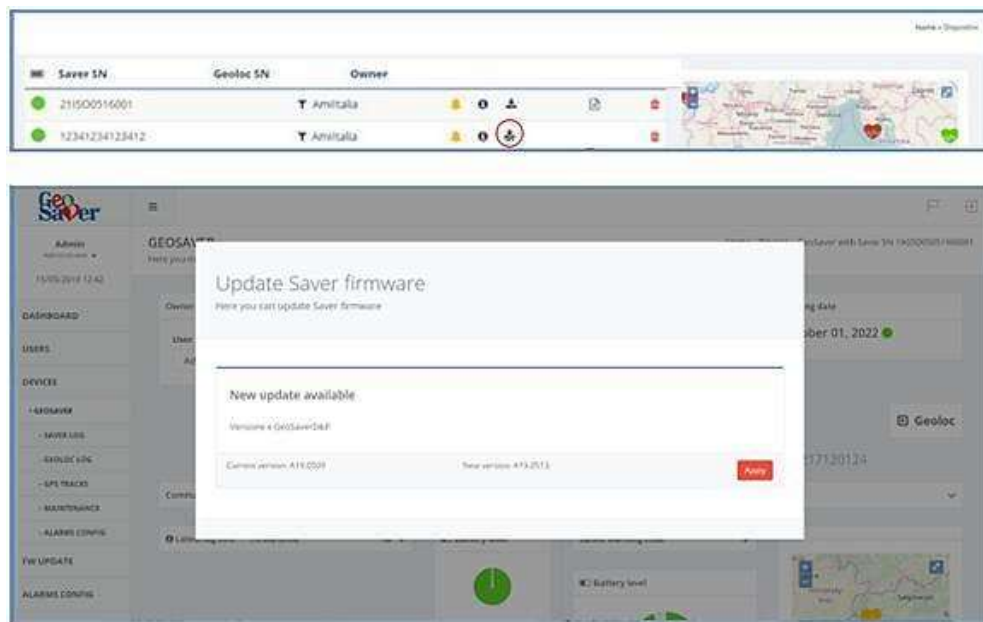
Press this key to save the set values

10.1.7 Firmware update

If AMI Italia makes new firmware available, the user is notified with an icon on the device table.


If a device has a less recent version of the firmware than the last one released on the line corresponding to the Geo Saver in question, an icon will be displayed, by clicking that you are taken to the device page. There the user can see which component needs an update: an icon the same as the one on the table is displayed next to the serial number of the device to be updated.

The update process begins with sending the command to the device: click on the aforementioned icon and follow the proposed steps.



10.2 Dashboard

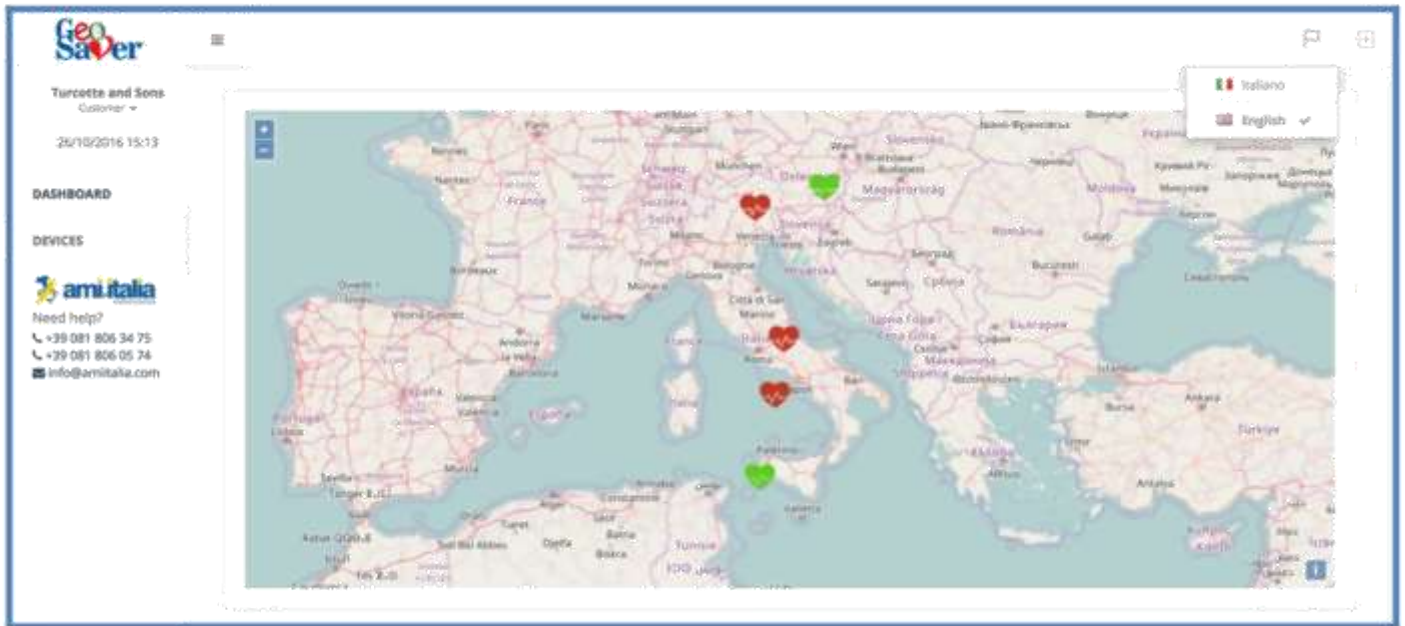After logging in to the cloud you are redirected to the Dashboard. In this page there is a map with all Geo Savers devices -linked to the current users- located through a coloured icon. The colour of the icon is representative of the state of the Geo Saver. Clicking one of these icons a small window appears with some of the main information related to that specific device.




10.3 Devices

In the section dedicated to the Devices, all the devices to which the current user has access are shown in a paged table. These devices can be filtered by:

- Membership user
- Range of serial numbers of the Saver
- Range of serial numbers of the Geoloc
- Geographic area

From the table it is possible to have quick access to some sections / actions by:

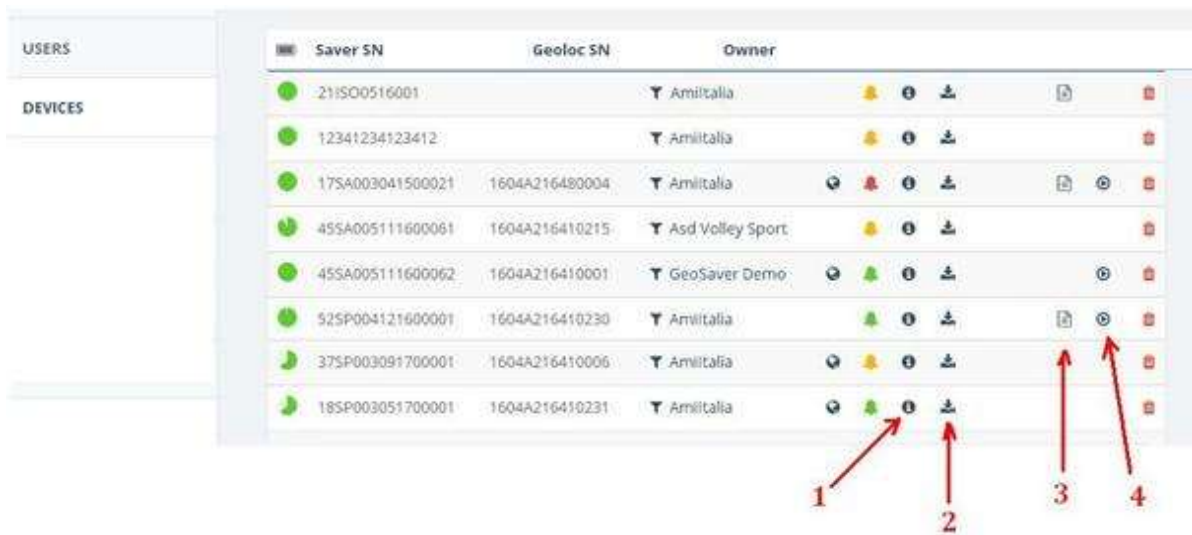1. Link to the page dedicated to the specific Geo Saver
2. Icon that notifies the presence of a new release (for at least one of the components - Saver, Geoloc, radio)
3. Link to the section that lists Complete Logs received from Geo Saver (where present)
4. Link to the section that lists the history of ECG streaming received from Geo Saver (where present)



USERS	Saver SN	Geoloc SN	Owner				
DEVICES	21500516001		Amitalia				
	12341234123412		Amitalia				
	17SA003041500021	1604A216480004	Amitalia				
	45SA005111600061	1604A216410215	Asd Volley Sport				
	45SA005111600062	1604A216410001	GeoSaver Demo				
	52SP004121600001	1604A216410230	Amitalia				
	37SP003091700001	1604A216410006	Amitalia				
	18SP003051700001	1604A216410231	Amitalia				

10.4 Geo Saver

The section dedicated to the specific Geo Saver can be divided into three parts

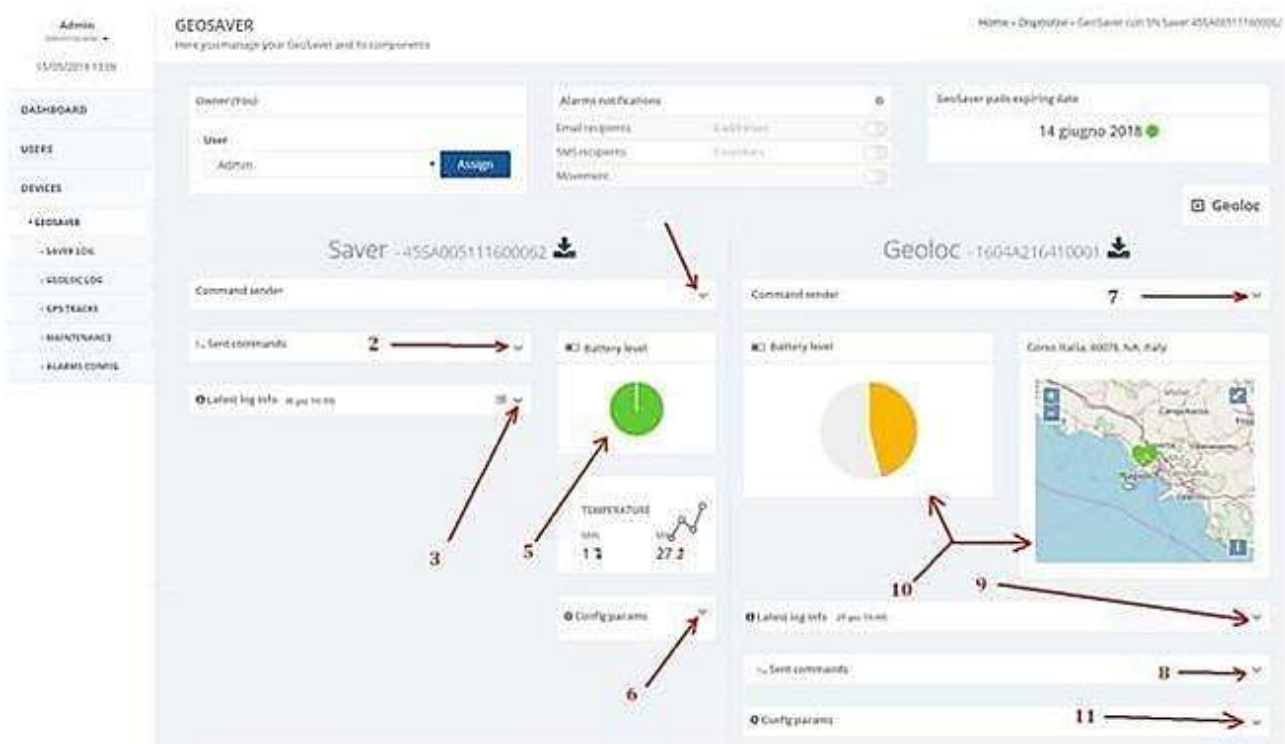
- At the top enabling / disabling of alarms
- On the left everything related to the Saver
- On the right everything related to the Geoloc (concealed)

In the area relating to the Saver we can distinguish:

1. Collapsible panel for sending a command
2. Collapsible panel with the last commands sent and relative completion status
3. Collapsible panel with the details of the last log received or Code of the last error
4. Code of the last warning (if present)
5. Pie chart with battery level
6. Collapsible panel with the last set of configurations

In the area relating to the Geoloc we can distinguish:


7. Collapsible panel for sending a command
8. Collapsible panel with the last commands sent and relative completion status
9. Collapsible panel with the details of the last log received or Code of the last error
10. Map showing the last position received or Pie chart with battery level
11. Collapsible panel with the last set of configurations



10.5 Saver Log

The section dedicated to logs received from the Saver presents a list of logs that can be expanded by clicking on them with the mouse.

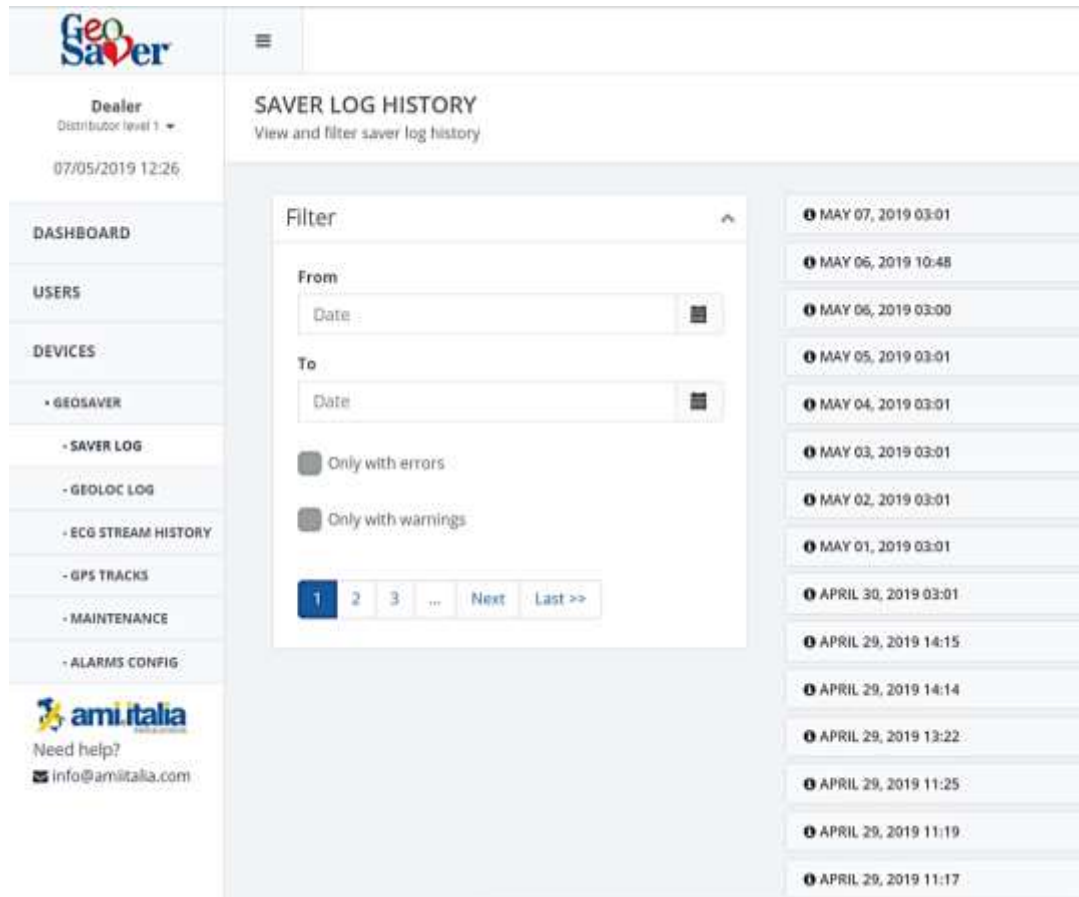On the left side of the page there is a box where it is possible to filter the list of logs arrived by arrival date. An additional filter can be added selecting only the logs with errors and / or warnings.



The screenshot shows the 'SAVER LOG HISTORY' interface. On the left is a navigation menu with options like 'DASHBOARD', 'USERS', 'DEVICES', and 'SAVER LOG'. The main area features a 'Filter' box with 'From' and 'To' date pickers, checkboxes for 'Only with errors' and 'Only with warnings', and a pagination control. To the right is a list of log entries with timestamps.

SAVER LOG HISTORY
View and filter saver log history

Filter

From
Date

To
Date

Only with errors

Only with warnings

1 2 3 ... Next Last >>

- MAY 07, 2019 03:01
- MAY 06, 2019 10:48
- MAY 06, 2019 03:00
- MAY 05, 2019 03:01
- MAY 04, 2019 03:01
- MAY 03, 2019 03:01
- MAY 02, 2019 03:01
- MAY 01, 2019 03:01
- APRIL 30, 2019 03:01
- APRIL 29, 2019 14:15
- APRIL 29, 2019 14:14
- APRIL 29, 2019 13:22
- APRIL 29, 2019 11:25
- APRIL 29, 2019 11:19
- APRIL 29, 2019 11:17

10.6 Geoloc Log

The section dedicated to logs received from the Geoloc presents a list of logs received which can be expanded by clicking on them with the mouse.

On the left side of the page there is a box where it is possible to filter the list of logs arrived by arrival date. An additional filter can be added by selecting only the logs with errors. There is also a map showing the position of the selected log or the last positions.

Dealer
Contributor level 1
07/05/2019 12:28

GEOLOG LOG HISTORY

View and filter geoloc log history

Filter

From
Date: [calendar icon]

To
Date: [calendar icon]

Only with errors

1 2 3 ... Next Last >>

amitalia
Need help?
info@amitalia.com

- MAY 07, 2019 03:01
- MAY 06, 2019 10:48
- MAY 06, 2019 03:00
- MAY 05, 2019 03:01
- MAY 04, 2019 03:01
- MAY 03, 2019 03:01
- MAY 02, 2019 03:01
- MAY 01, 2019 03:01
- APRIL 30, 2019 03:01
- APRIL 29, 2019 14:15
- APRIL 29, 2019 14:14
- APRIL 29, 2019 13:22
- APRIL 29, 2019 13:14
- APRIL 29, 2019 11:25
- APRIL 29, 2019 11:19

10.7 Streaming

During the streaming phase of a Geo Saver the link to this section will become visible, presenting:

- A graph showing the streaming in real time
- Some relevant data on the current stream
- The GPS position on a map

During the streaming, a Grawl notification is also presented with the link to this section and a red panel on the Geo Saver page.



10.8 List of saved ECG Streaming

The section dedicated to the list of ECG streams for the specific Geo Saver is made visible only in the presence of at least one file and presents:

- A graph to playback the streaming
- A filter from which to select the date of interest
- A table with the list of streaming resulting from application of the filter. From it you can start / stop playback or download the file.

ECG STREAM HISTORY Home » Devices » GeoSaver with Savitr SN 4350004101700051 » ECG Stream History
List of saver acquired ECG streams

Date

Filename	Acquired on	Duration	
2019-04-23_174331_4350004101700051.stream	April 23, 2019 17:46	00:03:35.396	
2019-03-26_180942_4350004101700051.stream	March 26, 2019 18:10	00:01:02.963	
2019-02-21_145848_4350004101700051.stream	February 21, 2019 14:58	00:01:03.388	
2018-11-09_122751_4350004101700051.stream	November 09, 2018 12:29	00:00:39.078	
2018-11-08_164536_4350004101700051.stream	November 08, 2018 16:46	00:00:43.853	

10.9 Tracks


The section dedicated to the saved tracks for the specific Geo Saver is made visible only in the presence of at least one file and presents:

- A map
- A filter to select the set of dates of interest
- A table with the list of tracks resulting from the application of the filter. By selecting one of these it is possible to view its layout on the map above.

TRACKS

View and filter device Tracks

Devices - GeoSaver with Saver SN 19G500501900001 - Tracks



Filter

From

Date:

To

Date:

Name	Filename	Acquired on		
2019-05-10_151004	1694A217120124_20190510_150356.gpx	May 10, 2019 15:03		
2019-05-10_145508	1694A217120124_20190510_140031.gpx	May 10, 2019 14:00		
2019-05-10_133831	1694A217120124_20190510_124501.gpx	May 10, 2019 12:45		
2019-05-09_173325	1694A217120124_20190509_174716.gpx	May 09, 2019 17:47		
2019-05-09_173408	1694A217120124_20190509_172738.gpx	May 09, 2019 17:27		

1
2
3
Next
Last >>

11.3 Cleaning

The structure of the *Geo Saver* defibrillator, including the connection port of the defibrillation electrodes, can be sanitized using a soft cloth dampened with one of the cleaning solutions listed below:

- a) Isopropyl alcohol (70% solution)
- b) Soap water
- c) Bleach (30 ml per litre of water)
- d) Detergents containing ammonia
- e) Detergents containing glutaraldehyde
- f) Oxygenated water



Do not immerse the *Geo Saver* in any liquid.

Do not use abrasive materials or detergents, strong solvents such as acetone or acetone-based detergents, and enzymatic detergents.

11.4 Preservation

The *Geo Saver* must be stored in a place where the environmental and safety conditions indicated in the table below are observed according to the temperature and humidity indicated in the chapter 12.2. If installed it is advisable to keep the device with the battery always inserted to allow it to carry out periodic self-diagnostic tests. For easy retrieval of the device in case of emergency, place it in easily accessible place and faced in a way that the control LEDs and the mini LCD are clearly visible.

	<p>Do not use, install or store the <i>Geo Saver</i> in conditions of temperature or humidity that exceed the ranges indicated in this user manual.</p>		<p>Do not install or store the <i>Geo Saver</i> in areas directly exposed to sunlight.</p>
	<p>Do not install or store the <i>Geo Saver</i> in areas subjected to sudden changes in temperature or humidity.</p>		<p>Do not install or store the <i>Geo Saver</i> near heat sources.</p>
	<p>Do not use, install or store the <i>Geo Saver</i> in places subjected to strong vibrations.</p>		<p>Do not use, install or store the <i>Geo Saver</i> in environments with high concentrations of flammable gases or anaesthetics.</p>
	<p>Do not install or store the <i>Geo Saver</i> in areas with a high concentration of dust.</p>		<p>The <i>Geo Saver</i> must be opened for maintenance only by A.M.I. Italia srl or by personnel authorized by the same.</p>

11.5 Troubleshooting Guide

The following table lists the symptoms, the possible causes and the possible corrective actions of the problems that may arise. For more information on the implementation of corrective actions, refer to the other sections of the user's manual. If the failure of the unit persists, request assistance.

Symptom	LED	Mini display	Possible cause	Corrective action
Device with battery installed does not switch on, the LED and the control display are both off.	OFF	OFF	The battery is totally discharged or faulty	Try replacing the battery. If the problem persists, call for assistance
			The device does not work	Ask for assistance
In standby the control LED flashes green but the mini display is off		OFF	The mini display is broken	Contact the assistance center
In standby the control LED is off but a "V" appears on the control mini-display.	OFF		The control LED is broken	Contact the assistance center
In standby the control LED flashes RED and a wrench appears on the control display.			During the daily self-test a critical error of the device was found	Contact a service center and report the error code.
In standby the control LED flashes GREEN / RED alternately and a wrench appears on the control display.	 		Very low battery Level <1% The device may turn off during use. <i>(see the relevant paragraph)</i>	Replace the battery
In the operating mode the voice command "Low battery" is issued.	 OFF		Low battery. 5% battery level. It is possible to use the device but the battery level is low <i>(see the relevant paragraph)</i>	Get a new battery and replace it as soon as possible.
During normal use the voice command "Battery low, Replace"	 		The battery is depleted. Level <1% The device may turn off during use. <i>(see the relevant paragraph)</i>	Avoid using the device if possible. Replace the battery
With the device turned on and after placing the PADS on the patient, the device continues to communicate: "Place Electrodes"	OFF		The PADS connector has not been inserted correctly or it has been removed	Insert the PADS connector in the appropriate compartment
			The PADS have been placed incorrectly	Correctly position the PADS on the patient's stripped chest. If necessary, remove the hair from the chest with a razor
			PADS are not working properly	Check the integrity and expiration of the PADS, replace them if necessary
The device turns on, the mini is on but no voice command is issued	OFF		The device's speaker does not work	Ask for assistance

12 Technical Specifications

The technical specifications of the *Geo Saver* defibrillator, its parts and accessories are shown below.

12.1 Physical characteristics

Category	Nominal specifications
Dimensions	29,5 x 23,0 x 11,5 cm
Weight	With battery Li-SOCl₂ (SAV-C1032): 2,54 Kg + Adult Pad (2,63 Kg)
	With battery Li-Ion (SAV-C1033): 2,57 Kg + Adult Pad (2,66 Kg)

12.2 Environmental requirements

Category	Nominal specifications
Temperature	Operational and standby: 0°C to 55°C (32°F a 131°F) Storage and transport: -40°C to 70°C (-40 a 158°F)
Relative humidity	Operational and standby: 10% to 95% (without condensation) Storage and transport: without humidity control (from -40°C to +5°C) up to 90% (from + 5 ° C to +35°C) with water vapor up to 50 hPa (from >35°C to +70°C)
Atmospheric pressure	Operating conditions: 620 hPa at 1060 hPa (altitude calculated min -382 mt and max 3955 mt)
Operating functional conditions	Normal use: Keep the AED device within the operating and standby ranges (not the storage and transport ranges) so that the device is ready for use. When starting from the inoperative conditions, let the device stabilize at the operating conditions for at least 2 hours, before the normal use.
Tolerance to impacts and falls	Complies with IEC/EN 60601-1 clause 21 (mechanical forces)
Sealing system	Complies with IEC/EN 60529 class IP56 standards; anti-spray, dustproof (with battery installed)
ESD (electrostatic shock)	Complies with IEC/EN 61000-4-2:2002 (3), Security level 4
EMC emissions / immunity	See chapter 13
Radio Equipment Directive	Directive 2014/53/UE

12.3 Reference regulations

Regulations and Directives	DIRECTIVE 2007/47/CE IEC/EN 60601-1 IEC/EN 60601-1-2 IEC/EN 60601-1-4 IEC/EN 60601-1-6 IEC/EN 60601-1-8 IEC/EN 60601-1-11 IEC/EN 60601-1-12 IEC/EN 60601-2-4 IEC/EN 60086-4 IEC/EN 60529 DIRECTIVE 2014/53/UE – RED ETSI EN 301 489-1, 7, 19, 52
----------------------------	--

12.4 Alarms Table

Priority	Cause	Visual signal
HIGH	Device ready to shock	LED Shock Icon flashing
HIGH	Low battery (<1% capacity)	Control LED flashing


12.5 Controls and indicators

Category	Nominal specifications
Buttons	ON / OFF button (device switching on and off) Shock button (to the defibrillation shock) “i” button INFORMATION
Visual Indicators	<ul style="list-style-type: none"> • Mini Display LCD control of device status • Device status control LED (bicolour RED / GREEN) • LED to place defibrillation PADS (2 LED Red) • LED do not touch patient (2 LED Red) • LED it is possible to touch patient (1 LED Red) • LED adult patient (1 LED Green) • LED paediatric patient (1 LED Green) • LED ON/OFF button (2 LED Green) • LED shock indicator AUTO (8 LED Red)
Sound Indicators	Multilingual voices for instructions during use of the device Acoustic signals of warnings and dangers
Speaker	Pre-set Volume (Emissions in compliance with IEC/EN 60601-2-4 point 6.1) Min. Variation 20% max 100% (60 dBA to 80 dBA ±3 dBA)
Microphone	Recording automatically activated when device is switched on
Streaming	ON / OFF setting from menu for sending ECG data to remote

12.6 Data memory

Category	Nominal specifications	
Internal Memory Capacity	6 hours of environmental audio recording, ECG tracing and events	
External memory (optional)	External SD / SDHC memory cards up to 8GB	
Archived data	AED1LOG.txt	Daily self-tests, Errors found, Device usage data, Device information
	AEDFILE.aed	Rescue events, voices and environmental noises, ECG tracing of rescue, Vital parameters of the patient analysed and detected by the Geo Saver
Data display	Via PC Saver View Express software (Microsoft Windows compatible)	

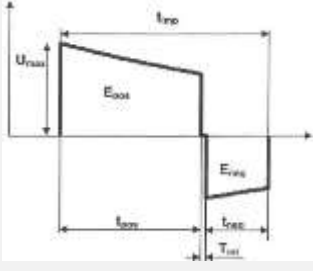12.7 Defibrillator

Category	Nominal specifications
Waveform 	<p>Biphasic Truncated Exponential (BTE)</p> <p>The waveform parameters are automatically adjusted according to the patient's impedance. In the graph on the left t_{pos} represents the duration of phase 1 (ms), t_{neg} represents the duration of phase 2 (ms), t_{int} is the delay between phases, U_{max} indicates the peak voltage, t_{imp} is the final voltage. In order to compensate for variations in the patient's impedance, the duration of each phase of the waveform is dynamically adjusted based on the charge delivered, as indicated in the paragraph following.</p>
Energy delivered (max) (Adults)	<p><i>Version 200J:</i> 200J nominal with a charge from 50 Ω</p> <p><i>Version 360J:</i> 350J nominal with a charge from 50 Ω</p>
Shock protocol (Adults)	<p><i>Version 200J:</i> Incremental: First: 150J – Subsequent: 200J</p> <p><i>Version 360J:</i> Incremental: First: 200J – Second: 250J – Subsequent: 350J</p>
Energy delivered (max) (Children)	<p><i>Version 200J:</i> 50 J nominal with a charge of 50 Ω</p> <p><i>Version 360J:</i> (when using defibrillation PADS SAV-C0016)</p>
Shock protocol (Children)	<p><i>Version 200J:</i> Fixed: First and Subsequent: 50J</p> <p><i>Version 360J:</i></p>
Charge control	Automatic through patient analysis system
Charge time (from the shock notice)	<p><i>Version 200J:</i> ≤ 9 SEC (according to IEC/EN 60601-2-4 §6.8.2 (7a)) (150J with new fully charged SAV-C1032 battery)</p> <p><i>Version 360J:</i> ≤ 15 SEC (according to IEC/EN 60601-2-4 §6.8.2 (7a)) (360J with new fully charged SAV-C1032 battery)</p>
Charge time (from the beginning of the analysis)	<p><i>Version 200J:</i> ≤ 15 SEC (according to IEC/EN 60601-2-4 §6.8.2 (8a)) (150J with new fully charged SAV-C1032 battery)</p> <p><i>Version 360J:</i> ≤ 21 SEC (according to IEC/EN 60601-2-4 §6.8.2 (8a)) (360J with new fully charged SAV-C1032 battery)</p>
Indication charge completed	<ul style="list-style-type: none"> • SHOCK button is flashing • Voice command “Press the red flashing button”
Shock delivery	The shock is delivered only through the SHOCK button
Disarmament	<p><i>Automatic:</i></p> <ul style="list-style-type: none"> • If the patient's analysis system considers the rhythm no longer shockable, or • If the operator has not pressed the SHOCK button within 15 seconds of completing charge, or • If defibrillation PADS have been removed from the patient or disconnected from the unit. <p><i>Manual:</i></p> <ul style="list-style-type: none"> • If the operator presses the OFF / DEACTIVATE button at any time to deactivate or switch off the appliance.
Shock detection vector	Through the defibrillation PADS (Lead II)
Isolation of the patient	Type BF

12.8 Efficiency of delivered energy


Impedance	Shock of 50 J (Paediatric)				Energy delivered (Joules)
	T _{pos} (ms)	T _{neg} (ms)	U _{max} (A)	Set Energy (J)	
25 Ohm	6,8	3,3	18,6	50	50,2
50 Ohm	7,2	3	12,3	50	49,2
75 Ohm	7,4	2,8	9,6	50	48,6
100 Ohm	7,5	2,7	8,1	50	48,4
125 Ohm	7,6	2,6	7,1	50	48,75
150 Ohm	7,7	2,5	6,4	50	48
175 Ohm	7,7	2,4	5,8	50	48,3

Impedance	Shock of 150 J				Energy delivered (Joules)
	T _{pos} (ms)	T _{neg} (ms)	U _{max} (A)	Set Energy (J)	
25 Ohm	4,6	5,6	43,8	150	147,2
50 Ohm	6,2	4	24,9	150	146,9
75 Ohm	6,8	3,3	18,4	150	147,15
100 Ohm	7,2	3	15	150	147,2
125 Ohm	7,4	2,8	13	150	146,5
150 Ohm	7,5	2,7	11,5	150	147
175 Ohm	7,6	2,6	10,4	150	147

Impedance	Shock of 200 J				Energy delivered (Joules)
	T _{pos} (ms)	T _{neg} (ms)	U _{max} (A)	Set Energy (J)	
25 Ohm	4,6	5,6	57,6	200	197,2
50 Ohm	6,1	4	28,8	200	196
75 Ohm	6,8	3,3	15,9	200	196,2
100 Ohm	7,2	3	17,3	200	196
125 Ohm	7,4	2,8	14,9	200	195,5
150 Ohm	7,5	2,7	13,2	200	195,3
175 Ohm	8,5	3	11,4	200	193,55

Impedance	Shock of 250 J				Energy delivered (Joules)
	T _{pos} (ms)	T _{neg} (ms)	U _{max} (A)	Set Energy (J)	
25 Ohm	4,6	5,6	56,6	250	246,4
50 Ohm	6,2	4	32,3	250	246,8
75 Ohm	6,8	3,3	23,7	250	244,95
100 Ohm	7,2	3	19,4	250	244,8
125 Ohm	8,4	3,4	15,8	250	241,75
150 Ohm	10	4	13,3	250	242,4
175 Ohm	11,5	4,6	11,4	250	241,15

Impedance	Shock of 350 J				Energy delivered (Joules)
	T _{pos} (ms)	T _{neg} (ms)	U _{max} (A)	Set Energy (J)	
25 Ohm	4,9	9,4	65,2	350	336
50 Ohm	7,2	6	36,6	350	341
75 Ohm	9,5	6,9	25,4	350	339,3
100 Ohm	12	8,2	19,4	350	339
125 Ohm	14,4	9,5	15,8	350	338,5
150 Ohm	16,9	10,9	13,3	350	339
175 Ohm	18,9	11,5	11,4	350	331,8

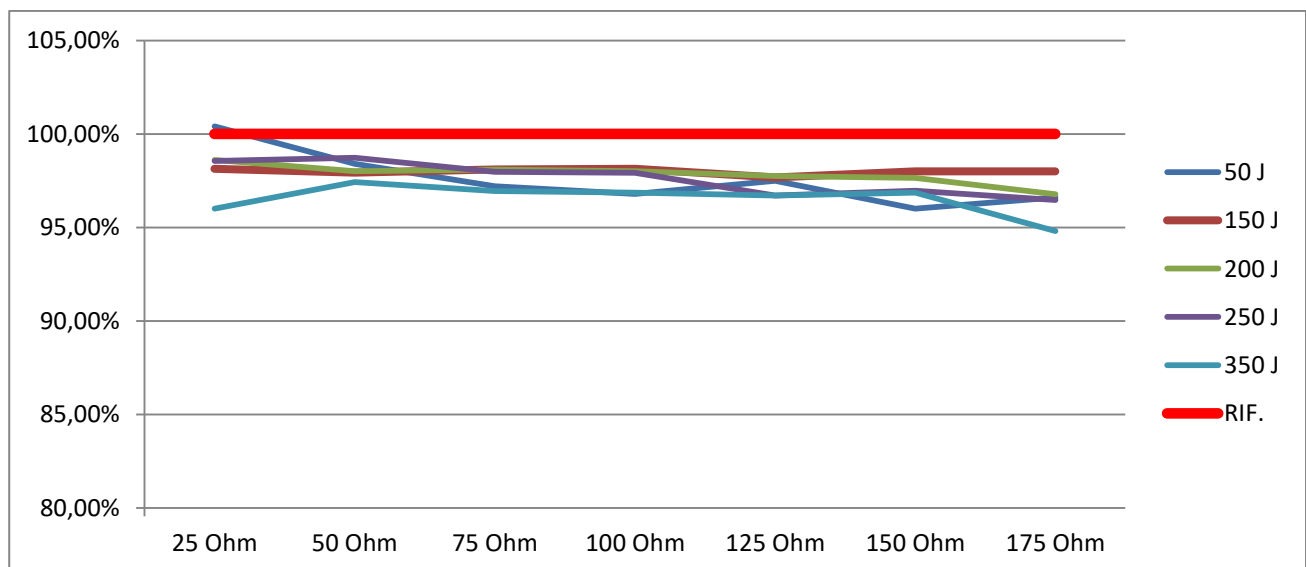

Efficiency of the energy supplied graph

12.9 Patient analysis system

Category	Nominal specifications
Function	Determines the patient's impedance and evaluates the ECG rhythm and signal quality to determine whether or not the shock delivery is appropriate.
Impedance range	20- 200 Ω
ECG analysis time	≥ 4 seconds (with new fully charged battery) respecting IEC/EN 60601-2-4
Sensibility	97% Respects the guidelines IEC/EN 60601-2-4 2002(3) sources AHADB, MITDB and EDB
Specificity	99% Respects the guidelines IEC/EN 60601-2-4 2002(3) sources AHADB, MITDB and EDB
Shockable rhythms	<p>If used on a patient who has the characteristics listed in the usage criteria, the <i>Geo Saver</i> defibrillator is designed to recommend a defibrillating shock when it detects the right impedance and when the following situations occur:</p> <p>Ventricular Fibrillation peak-to-peak amplitude at least 200μVolts Ventricular Tachycardia with cardiac rhythm frequency min. 180 bpm and peak-to-peak amplitude at least 200μVolts (including ventricular flutter and polymorphic ventricular tachycardia)</p>
Non-shockable rhythms	The <i>Geo Saver</i> is designed to not recommend shock with all other rhythms, including: normal sinusoidal rhythm, moderate ventricular fibrillation (<200 μ Volts), some slow ventricular tachycardia and asystoles.

12.10 ECG Analysis Function

ECG rhythm	Dimension Test sample	Objective	Detected value
Shockable rhythm Ventricular Fibrillation (VF)	500	Sensibility > 90%	98%
Shockable rhythm Ventricular Tachycardia (VT, bpm>140)	600	Sensibility > 75%	92%
Non-shockable rhythm Normal sinus rhythm	1500	Specificity > 99%	100%
Non-shockable rhythm Asystole	30	Specificity > 95%	100%
Untreatable rhythm generic AF, SVT, PVC	30	Specificity > 95%	100%
Positive predictive values			97.1%
False positives			4.1%

12.11 Non rechargeable battery

Category	Nominal specifications						
REF (Model)	SAV-C1032						
Type	Li-SOCl ₂ (lithium-thionyl chloride) disposable, non-rechargeable						
Voltage	28,8 VDC – 3500 mAh						
Microphone	Active if memory card is installed to record voices and environmental noises						
Capacity	<table border="0"> <tr> <td>Standard 200J</td> <td>300 cycles of complete rescues (shocks at 200J. and CPR)*</td> </tr> <tr> <td>Power 360J</td> <td>200 cycles of complete rescues (shocks at 360J. and CPR)*</td> </tr> <tr> <td>ECG analysis</td> <td>35 continuous hours*</td> </tr> </table>	Standard 200J	300 cycles of complete rescues (shocks at 200J. and CPR)*	Power 360J	200 cycles of complete rescues (shocks at 360J. and CPR)*	ECG analysis	35 continuous hours*
Standard 200J	300 cycles of complete rescues (shocks at 200J. and CPR)*						
Power 360J	200 cycles of complete rescues (shocks at 360J. and CPR)*						
ECG analysis	35 continuous hours*						
Duration in Standby (installed battery)	4 years if installed in the AED, assuming an activation test, daily self-tests without turning on the AED						

* New and fully charged battery at a constant temperature of 20°C and relative humidity without condensation 45%

12.12 Rechargeable battery

Category	Nominal specifications	
REF (Model)	SAV-C1033	
Type	Li ion (lithium ions) Rechargeable	
Voltage	21,6 VDC - 2100 mAh	
Capacity	<i>Standard 200J</i>	250 continuous shock*
	<i>Power 360J</i>	160 continuous shock *
	<i>ECG analysis</i>	21 continuous hours*
Charging time	≤ 2,5 hours with charging station SAV-C1035*	
Shelf Life	2 years or 300 charge/shock cycles (<i>the one that occurs first</i>) *	

* New and fully charged battery at a constant temperature of 20°C and relative humidity without condensation 45%

12.13 Internal back-up battery

Category	Nominal specifications
Type	Battery Coin Cell (LiMnO ₂)
Purpose	Maintaining configuration data (date / time, etc)
Voltage	3 VDC
Duration	Maintains data for 3 years (without external battery) Maintains data for 6 years (with external battery inserted within 12 months)

12.14 Rechargeable battery charger

Category	Nominal specifications	
REF (Model)	SAV-C1035	
Charge control	LED multicolour red green (see paragraph 5.1.6)	
Power supply	<i>Input</i>	15Vdc-2.67A / 12Vdc-5.5A
	<i>Output</i>	26VDC – 1,5A
	<i>Absorption</i>	40W/66W
AC/DC Adapter	<i>Model</i>	MeanWell GS40A15-P1J
	<i>Identification code</i>	SAV-C1037
	<i>Input</i>	100-240VAC – 50/60Hz – 1.5A
	<i>Output</i>	15V – 2.67A
	<i>Absorption</i>	40W

12.15 Defibrillation PADS

Category	ADULTS	CHILDREN
REF (Model)	SAV-C0846	SAV-C0016
Series	Cable and connector outside the envelope	Cable PADS and connector inside the envelope
Patient range	Adults age >8 years or weight > 25Kg	Children age 1 - 8 years or weight < 25Kg
Intended use	Disposable	
No. of shock tolerated	50 shocks at 360J	
Support material	Medical FOAM, thickness 1 mm	
Conductive gel	Low impedance conductive adhesive gel	
Total area (<i>for pad</i>)	136 cm ²	75 cm ²
Active area (<i>for pad</i>)	94 cm ²	40 cm ²
Conductive material	Metal foil	
Connection	Anti-shock safety connector	
Cable length	120 cm (normally)	

12.16 Timing of Shock cycles

Charging time performance in accordance with 60601-2-4 (201.101)	Specific	Result
The maximum time between the beginning of the ECG rhythm analysis and the completion of the charge at maximum energy	< 30 seconds	OK
The maximum time from turn on to completion of the charge at maximum energy	< 40 seconds	OK

12.17 Geoloc Module

Geoloc Module	
Frequency	GSM: 850, 900, 1800, 1900 MHz; UMTS: 900, 2100 MHz GPS: 1575, 1600 MHz
Performance	<ul style="list-style-type: none"> • Geo-location (geographical position of the device in real time) • Remote control (tracking, anti-theft and device configuration) • Telemetry (real-time acquisition of device status, parameters and alarms) • Remote assistance (real-time phone call and ECG data streaming)

12.18 Geoloc Module Battery Type

Geoloc Module Battery Type	
Type	Li-SOCl ₂ Disposable code SAV-C1038
Autonomy	4 years after installation on AED with a battery insertion test and daily self-test but without activating AED*
Shelf-Life	5 years if stored in its original packaging*
Type	Li-Ion Accumulator (rechargeable), code SAV-C1039
Recharging time	2 hours with the charging station code SAV-C1040 *
Autonomy	2 years after installation on AED and running only the daily self-diagnosis, but without activating AED*
Shelf-Life	2 years or 300 charge cycles*

*New and fully charged battery at constant temperature of 20°C and non-condensing relative humidity 45%.

13 Compliance with electromagnetic emission standards

The following paragraphs will specify the compliance with electromagnetic emission standards:

- Guidelines and manufacturer's declaration - Electromagnetic emissions
- Guidelines and manufacturer's declaration - Electromagnetic immunity
- Recommended distances between portable and mobile radiofrequency communication equipment and the AED

13.1 Guidelines and manufacturer's declaration - Electromagnetic emissions

The *Geo Saver* was designed to be used in electromagnetic environments with features listed below. The customer or the user of the *Geo Saver* must ensure that it is used in such an environment.

Emissions test	Compliance	Electromagnetic environment - Guidelines
RF Emissions CISPR 11	Group 1	The AED uses RF energy only for its internal operation. Therefore its RF emissions are very low and are unlikely to interfere with nearby electronic equipment.
RF Emissions CISPR 11	Class B	The AED can be used in any building, including those for residential use and those directly connected to the public low-voltage power supply network that supplies residential buildings.
Harmonic Emissions IEC/EN 61000-3-2	Not applicable	
Voltage fluctuations / flicker IEC/EN 61000-3-3	Not applicable	

13.2 Guidelines and manufacturer's declaration - Electromagnetic immunity

The *Geo Saver* was designed to be used in electromagnetic environments with features listed below. The customer or the user of the *Geo Saver* must ensure that it is used in such an environment.

Immunity test	Test level IEC 60601-1	Compliance level	Electromagnetic environment - Guidelines
Electrostatic shock (ESD) IEC/EN 61000-4-2	±6 kV contact	±6 kV contact	Floors must be wood, concrete or ceramic tiles. If the floors are covered with synthetic material, the relative humidity must be at least 30%.
	±8 kV air	±8 kV air	
Fast transients / bursts IEC/EN 61000-4-4	±2 kV by electricity networks	Not applicable	
	±1 kV by input / output networks	±1 kV by input / output networks	
IEC/EN 61000-4-11	< 5% U_T (> 95% dip in U_T) for 0,5 cycles 40% U_T (60% dip in U_T) for 5 cycles 70% U_T (30% dip in U_T) for 25 cycles < 5% U_T (>95% dip in U_T) for 5 seconds	Not applicable	

Immunity test	Test level IEC/EN 60601-1	Compliance level	Electromagnetic environment - Guidelines
Supply frequency (magnetic field) 50/60 Hz IEC/EN 61000-4-8	3 A/m	80 A/m	Power frequency magnetic fields must be at levels no higher than those of stations located in typical heavy industrial applications, power plants and control rooms of high voltage substations.
Note: U_T is the main AC current before the test level is applied			
RF conducted	3 Vrms	Not applicable	
IEC/EN 61000-4-6	from 150 kHz to 80 MHz outside of ranges ISM ^a 10 Vrms from 150 kHz to 80 MHz inside the ranges ISM ^a	Not applicable	
RF radiated IEC/EN 61000-4-3	10 V/m from 80 MHz to 2,5 GHz	10 V/m	<p>The distance between portable and mobile RF communications equipment in use and any part of the AED, including cables, must never be less than the recommended separation distance calculated based on the equation applicable to the transmitter frequency.</p> <p>Recommended separation distance</p> $d = 1.2\sqrt{P}$ from 80 MHz to 800 MHz $d = 2.3\sqrt{P}$ from 800 MHz to 2,5 GHz Where P is the transmitter's maximum output power range in watts (W) according to the transmitter manufacturer's data and d is the recommended distance in meters (m) ^b . The field strengths of fixed radiofrequency transmitters, as determined by an investigation in electromagnetic sites, ^c should be less than the compliance level in each frequency range. ^d Interference may occur near the devices marked with this symbol. 
NOTE 1	From 80 MHz to 800 MHz, the higher frequency range applies.		
NOTE 1	These guidelines may not apply in all situations. Electromagnetic propagation is influenced by absorption and reflection from structures, objects and people.		
a	The ISM bands (industrial, scientific and medical) between 150 kHz and 80 MHz are from 6.765 MHz to 6.795 MHz; from 13.553 MHz to 13.567 MHz; from 26.957 MHz to 27.283 MHz; and from 40.66 to 40.70 MHz.		
b	The levels of compliance in the ISM bands between 150 kHz and 80 MHz and between 80 MHz and 2.5 GHz are designed to reduce the possibility of interference in the event that portable and mobile communication devices are inadvertently approached to the area where the patient is found. For this reason, an additional factor of 10/3 is added to the calculation of the recommended separation distance for transmitters whose frequencies fall within these ranges.		
c	It is not possible to accurately predict theoretically the field strengths of fixed transmitters, such as base stations for radio (cellular / cordless) telephones and mobile radios, amateur radio, AM and FM radio and TV. To assess the electromagnetic environment with fixed RF transmitters, consider conducting an electromagnetic site survey. If the power of the fields measured in the location in which the AED is used exceeds the specific RF compliance level mentioned above, it will be necessary to monitor the AED to verify its correct functioning. If operating anomalies are observed, it may be necessary to take corrective measures, for example by moving or reorienting the AED.		
d	Over the frequency range between 150 kHz and 80 MHz, field strengths must be less than 1 V / m.		

13.3 Recommended separation distance between portable and mobile RF communication equipment and Geo Saver device

The **Geo Saver** must be used in an electromagnetic environment in which radiated RF interference is controlled. The customer or the operator of the **Geo Saver** can help prevent electromagnetic interference by maintaining the minimum distances recommended below, between the portable and mobile RF communications equipment (transmitters) and the **Geo Saver**, based on the maximum output power of the devices of communication.

Maximum transmitter power output rate W	Separation distance according to the transmitter frequency m			
	From 150kHz to 80 MHz outside the ISM bands	From 150kHz to 80 MHz inside the ISM bands	From 80 MHz to 800 MHz	From 800 MHz to 2,5 GHz
	$d = 1.2\sqrt{P}$	$d = 1.2\sqrt{P}$	$d = 1.2\sqrt{P}$	$d = 2.3\sqrt{P}$
0.01	0,12 m	0,12 m	0,12 m	0,23 m
0.1	0,37 m	0,38 m	0,38 m	0,73 m
1	1,12 m	1,2 m	1,2 m	2,3 m
10	3,7 m	3,8 m	3,8 m	7,3 m
100	12 m	12 m	12 m	23 m
For transmitters rated at a maximum power not listed above, the separation distance "d" in meters (m) can be determined using the equation applicable to the transmitter frequency, where P represents the maximum power produced by the watt transmitter (W) according to the transmitter manufacturer.				
NOTE 1:	At 80 MHz and 800 MHz, the separation distance applied is that used for high frequency ranges.			
NOTE 2:	The ISM frequency bands (for industrial, scientific and medical applications) between 150 kHz and 80 MHz are 6,765 MHz up to 6,795 MHz; 13,553 MHz up to 13,567 MHz; 26,957 MHz up to 27,283 MHz and 40,66 MHz up to 40,70 MHz			
NOTE 3:	An additional factor of 10/3 is used in the calculation of the recommended separation distance for transmitters in the ISM frequency bands between 150 kHz and 80 MHz and in the frequency range from 80 MHz to 2.5 GHz to decrease the possibility that a Mobile / portable equipment may interfere if inadvertently brought into the patient's area.			
NOTE 4:	These guidelines may not be applicable in all situations. Electromagnetic diffusion is influenced by the absorption and reflection of structures, objects and people.			

14 Simbology

	ILCOR Universal Symbols for AED
	High Voltage Electrical Hazard
	General Notices: Refer to the consultation of accompanying documents before using the appliance
	Type BF, Defibrillation-proof Equipment
	Do not expose to high temperatures or flames
	Do not recharge
	Do not open
	Do not destroy or damage it
	Do not use it in puddles of water
	Read the User Manual
	Battery recycling
	Follow local waste regulations
	Fragile
	Store in a dry place
	Do not expose to direct sunlight
	Risk of electric shock do not open
	Type CF applied part

	IMQ brand
	CE mark with identification number
	Degree of protection of the appliance against dust and water (including battery)
	Serial Number
	Manufacturing date
	Lot Number (LOT)
	Expiration date
	Model identifier
	Manufacturer Name
	Absence of latex
	Single use, do not reuse
	Not Sterile
	External directions on the box
	This side up
	Temperature Limits
	Only stack up to 6 cartons in height
	Non-ionizing electromagnetic radiation

15 Certifications

15.1 EC Certificate

Mod. 4606/0



CERTIFICATO CE

Certificato n. 1104/MDD

Dichiarazione di approvazione del sistema qualità

(Sistema completo di garanzia qualità)

Visto l'esito delle verifiche condotte in conformità all'Allegato II, con l'esclusione del punto 4, della direttiva 93/42/CEE e s.m.i., si dichiara che la ditta:

A.M.I. ITALIA S.R.L.

80143 NAPOLI (NA) - VIA G. PORZIO CENTRO DIREZIONALE IS.G2 (ITA) - Italy

mantiene nello stabilimento di:

A.M.I. INTERNATIONAL KFT - 2000 SZENTENDRE - KOZUZO u. 5/A (HUN) - Hungary

80010 QUARTO (NA) - VIA CUPA REGINELLA 15A (ITA) - Italy

un sistema qualità che assicura la conformità dei seguenti prodotti:

Defibrillatore cardiaco esterno

Modd. come da documento "Defibrillatore Cardiaco Esterno" Rev.0 del 09/11/2018; valido solo se provvisto del timbro IMQ.

ai requisiti essenziali della direttiva suddetta ad essi applicabili (in tutte le fasi dalla progettazione al controllo finale) ed è sottoposta alla sorveglianza prevista dal punto 5 dell'Allegato II. Per i dispositivi in classe III questo certificato è valido solamente con il relativo certificato di esame CE della progettazione di Allegato II.4.

Riferimento pratiche IMQ:
 10AJ00006; 10AJ00117; COMEDCONMHDM110027747-01; 10EN00018; 10AC00009; DM17-0009799-01;
 DM17-0018806; DM17-0020656-01; DM18-0023720-01; DM18-0032037-01; DM19-0034531-01.

Questa Dichiarazione di approvazione è rilasciata dall'IMQ S.p.A. quale organismo notificato per la direttiva 93/42/CEE e s.m.i. Il numero identificativo dell'IMQ S.p.A. quale organismo notificato è: 0051.

Emesso il:	2008-02-18
Data aggiornamento:	2019-02-22
Sostituisce:	2018-11-15
Data scadenza:	2023-02-15




Questa Dichiarazione di approvazione è soggetta alle condizioni previste dall'IMQ nel "Regolamento per la certificazione CE dei dispositivi medici - Marcatura CE - Direttiva 93/42/CEE".

IMQ S.p.A. | I-20138 Milano
 | Via Quintiliano 43 |
 www.imq.it



EC CERTIFICATE

Certificate No 1104/MDD

Full Quality Assurance System Approval Certificate

On the basis of our examination carried out according to Annex II, excluding section 4, of the Directive 93/42/EEC and its revised version, we hereby certify that:

A.M.I. ITALIA S.R.L.

80143 NAPOLI (NA) - VIA G. PORZIO CENTRO DIREZIONALE IS.G2 (ITA) - Italy

manages in the factory of:

A.M.I. INTERNATIONAL KFT - 2000 SZENTENDRE - KOZUZO u. 5/A (HUN) - Hungary

80010 QUARTO (NA) - VIA CUPA REGINELLA 15A (ITA) - Italy

a quality assurance system ensuring the conformity of the following products:

External cardiac defibrillator

Type ref. as to Document "Defibrillatore Cardiaco Esterno" Rev.0 dated 2018/11/09; valid only if provided with IMQ mark.

with the relevant essential requirements of the aforementioned directive (from design to final inspection and testing) and it is subject to surveillance as specified in section 5 of Annex II. For class III devices, this certificate is valid only with the relevant EC Design-Examination Certificate of Annex II.4.

Reference to IMQ files Nos:

10AJ00006; 10AJ00117; COMEDCONMHDM110027747-01; 10EN00018; 10AO00009; DM17-0009799-01; DM17-0018806; DM17-0020656-01; DM18-0023720-01; DM18-0032037-01; DM19-0034531-01.

This Approval Certificate is issued by IMQ S.p.A. as Notified Body for the Directive 93/42/EEC and its revised version. Notified Body notified to European Commission under number: 0051.

Date: 2008-02-18
 Updated: 2019-02-22
 Substitution Date: 2018-11-15
 Expiry Date: 2023-02-15

This Approval Certificate is subjected to the provisions laid down in the "IMQ regulation for the certification of Medical Devices - CE Marking - Directive 93/42/EEC".

This is a translation of the Italian text, which prevails in case of doubts

IMQ S.p.A. | I-20138 Milano
 | Via Quintiliano 43 |
 www.imq.it

15.2 IMQ Brand



IMQ S.p.A. - Società con Socio Unico
I-20138 Milano - via Quintilano, 43
tel. 0250731 (p.a.) - fax 0250991500
e-mail: info@imq.it - www.imq.it

Rea Milano 1595884
Registro Imprese Milano 12898410159
C.F./P.I. 12898410159
Capitale Sociale € 4.000.000

CA10.00185
SN.I000XN

PID:
10010024
CID:
CN.I0005Y

Certificato di approvazione

Approval certificate



IMQ, ente di certificazione accreditato, autorizza la ditta

IMQ, accredited certification body, grants to

PRD N° 005B

Membro degli Accordi di Mutuo Riconoscimento EA, IAF e ILAC
Signatory of EA, IAF and ILAC Mutual Recognition Agreements

A.M.I. ITALIA S.R.L.
VIA G. PORZIO CENTRO DIREZIONALE IS.G2
80143 NAPOLI NA
IT - Italy

all'uso del marchio

the licence to use the mark

IMQ

Il presente certificato è soggetto alle condizioni previste nel Regolamento "MARCHI IMQ - Regolamento per la certificazione di prodotti" ed è relativo ai prodotti descritti nell'Allegato al presente certificato.



per i seguenti prodotti

for the following products

Defibrillatori cardiaci esterni
(Modd.: SAVER ONE; SAVER ONE D; SAVER ONE P; GEO SAVER; GEO SAVER D; GEO SAVER P)

External cardiac defibrillators
(Models: SAVER ONE; SAVER ONE D; SAVER ONE P; GEO SAVER; GEO SAVER D; GEO SAVER P)

This certificate is subjected to the conditions foreseen by Rules "IMQ MARKS - RULES for product certification" and is relevant to the products listed in the annex to this certificate.

Emissa il / Issued on **2008-09-25**
Aggiornato il / Updated on **2019-03-04**
Sostituisce / Replaces **2014-03-18**

Stefano D. M.
IMQ S.p.A.

16 Geo Saver Series Defibrillator Warranty

1 Warranty Restriction

A.M.I. Italia S.r.l guarantees the original purchasers that its Geo Saver series defibrillators and related accessories and batteries are free from any material or manufacturing defect according to the terms and conditions of this restrictive warranty. The original purchaser is considered to be the final user of the product purchased. This limited warranty is granted only to the original purchaser of the Geo Saver series defibrillator of A.M.I. Italia S.r.l and is not transferable or assignable to third parties.

The Geo Saver One Series defibrillators are as follows:

Geo Saver Semi-Automatic (code SGS-B0988 /SGS-B0989)

Geo Saver Automatic (code SGA-B0990/ SGA-B0991)

Geo Saver D (code SGD-B0992 / SGD-B0993)

Geo Saver P (code SGP-B0994 / SGP-B0995)

2 Duration

A.M.I. Italia S.r.l guarantees the original purchaser of Geo Saver series defibrillators, starting from the date of dispatch* of the warranty validation form (to A.M.I. Italia S.r.l) or starting from 30 (thirty) days from the date of shipment from A.M.I. Italia S.r.l, the one that occurs chronologically first; defibrillators have a typical life expectancy of about 10 years. The guarantee offered by A.M.I. Italia S.r.l covers a period equal to:

- **AED Geo Saver Series** have a six (6) year warranty
- **Non-rechargeable batteries Li-SOCl₂** (SAV-C1032) if installed in the AED and in Standby mode they are guaranteed for 4 (four) years assuming a battery activation test, daily self-tests, without the AED being switched on at the following environmental conditions temperature (20 ° C) and humidity S / C (45 %)
- **Rechargeable batteries Li-Ion** (SAV-C1033) are guaranteed for two (2) years from the date of production only if the temperature conditions (temperature 20 ° C) and humidity (45%) are met and if they are recharged at least one (1) time every four (4) months
- **The disposable pads** guaranteed until their expiration date.
- All **other accessories** are guaranteed for six (6) months starting 30 days after the original shipping date from our warehouse.

**The date shown on the registered letter with return receipt will still be valid*

3 Procedure

Please complete (in its entirety) the limited warranty validation form and send it by post (Registered letter A / R) to A.M.I. Italia S.r.l. The date shown on the A / R recommendation will prevail. You will find the Warranty validation form attached to the user manual or inside the original packaging of the Geo Saver series defibrillator. In the event that a defect covered by this warranty is found, the original purchaser must contact the reference retailer or an authorized A.M.I. Italia S.r.l.

A.M.I. Italia Srl reserves the exclusive right to repair or replace the product.

4 Exclusions

This warranty does not cover non-conformities subsequent to purchase, such as those caused by accidents, modifications, negligence, incorrect use or abuse, non-compliance with procedures or hazards, or warnings or cautions described in the user manual, failure to perform a reasonable and adequate maintenance, incorrect installation, replacement of parts and accessories that do not comply with the specifications provided by AMI Italia Srl, any modifications made to the device and in general all subsequent non-conformities deriving from failure to comply with the provisions contained in the user manual.

This warranty does not cover, as it does not constitute cases of original non-conformity, the normal wear and tear of components subject to decay during use such as buttons, LEDs and battery contacts. This warranty will also be automatically invalidated in one of the following cases:

- if the Geo Saver AED series serial number is modified, deleted, rendered illegible or otherwise tampered with;
- if the warranty seal (opening of the device) on the Geo Saver DAE is removed;
- in case the commercial name of the product or manufacturer is covered, modified or cancelled

Finally, this warranty does not apply to used Geo Saver series AEDs sold, in which case the warranty must be offered by the reseller of the used product with the exclusion of any liability, even indirect, borne by A.M.I. Italia S.r.l

5 Damage

Except as explicitly provided in this warranty, A.M.I. Italia S.r.l, WILL NOT BE LIABLE FOR ANY INCIDENTAL OR INDIRECT DAMAGES ARISING FROM THE USE OF THE GEO SAVER SERIES DEFIBRILLATOR OR CLAIMS IN VIRTU OF THIS AGREEMENT, WETHER THE CLAIM REFERS TO THIS AGREEMENT, TO ILLEGAL OR OTHERWISE. The warranty statements mentioned are exclusive and replace any other remedy. Some states do not allow the exclusion or limitation of incidental or indirect damages, so the above limitation or exclusion may not be relevant.

6 Waiver

ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE AND ALL IMPLIED WARRANTIES ARISING OUT OF NEGOTIATIONS, USE OR BUSINESS CONSUMPTIONS, BY STATUTE OR OTHERWISE ARE STRICTLY LIMITED TO THE TERMS OF THIS WRITTEN WARRANTY. This warranty will be your sole and exclusive buyer's remedy for this purchase. In the event of an alleged violation of any guarantee or legal action brought by the original purchaser for alleged negligence or other unlawful conduct by A.M.I. Italia S.r.l, the sole and exclusive remedy of the original purchaser will be constituted by the repair or replacement of the resulting defective materials, based on what was previously established. No retailer or agent or employee of A.M.I. Italia S.r.l is authorized to make changes, extensions or additions to this warranty.

7 Territorial limits

This warranty is valid for products purchased in one of the countries of the European Union or in the countries in which the EU laws and regulations apply.

8 Warning

Install, use and maintain Geo Saver series defibrillators by A.M.I. Italia S.r.l in absolute compliance with the instructions contained in the user manual

9 Other rights

This limited warranty guarantees the original purchaser specific legal rights; any other rights may vary depending on the state of belonging.

10 Applicable laws

Any dispute relating to this agreement or arising from the use of Geo Saver series defibrillators by A.M.I. Italia S.r.l will be governed by Italian law, at the Court of Naples, Italy

17 Product registration

In order to guarantee a correct and rapid traceability of the product sold, we ask you to complete the form below and send it by fax or registered letter to A.M.I. Italia S.r.l.

Warranty Card
Scheda di Garanzia
Carte de Garantie
Tarjeta de Garantía

GeoSaver Series

Complete AEDs range to save lives

<p>Device Model Modello del Dispositivo Modèle du Dispositif Modelo de Dispositivo</p>	<p>Serial Number (see label on the back) Numero di Serie (vedere l'etichetta sul retro) Numéro de Série (voir l'étiquette sur le dos) Número de Serie (ver la etiqueta en la parte posterior)</p>	<p>Date of Purchase Data d'Acquisto Date d'Achat Fecha de la Compra</p>	
<p>End User's Name Nome dell'Utente Finale Nom de l'Utilisateur Final Nombre del Usuario Final</p>	<p>Address Indirizzo Adresse Dirección</p>	<p>Country Paese Pays País</p>	
<p>City Città Ville Ciudad</p>	<p>State/Province/Region Stato/Provincia/Regione Etat/Province/Région Estado/Provincia/Región</p>	<p>Postal Code Codice Postale Code Postal Código Postal</p>	
<p>Telephone Number Numero di Telefono Numéro de Téléphone Número Telefónico</p>	<p>Fax Number Numero di Fax Numéro de Fax Número de Fax</p>	<p>Email Address Indirizzo Email Adresse Email Correo Electrónico</p>	
<p>Seller's Company Name Ragione Sociale del Venditore Nom du Vendeur Nombre del Vendedor</p>	<p>Seller's Country Paese del Venditore Pays du Vendeur País del Vendedor</p>	<p>Seller's Phone Number Telefono del Venditore Téléphone du Vendeur Teléfono del Vendedor</p>	



GeoSaver®

AED_S

